



Steven M. Ruben
Appl. No. 10/662,429

Account Name: Steven M. Ruben
Lab Notebook

Department MOL. BIO.
Subject 25796 - 4/8/96
Name ANN KIM # 12
Address _____
National Brand
Computation Notebook
11 3/4" x 9 1/4", 4 x 4 Quad., 75 Sheets 43-648

0 73333 43648 8
 **AVERY DENNISON**
Office Products
Chicopee, MA 01022

BEST AVAILABLE COPY

Ruben EXHIBIT #95

2

Department MOL. BIO.
Subject 2/5/96 - 4/8/96
Name ANN KIM #12
Address _____

National Brand
Computation Notebook
11 3/4" x 9 1/4", 4 x 4 Quad., 75 Sheets **43-648**



0 73333 43648 8

 **EVERY DENNISON**
Office Products
Chicopee, MA 01022

Ruben EXHIBIT 2095
Ruben v. Wiley et al.
Interference No. 105,077
RX 2095

2/6/96

DNA	4	Inoculate 37°C O/N
10x#2	5	
H ₂ O	40	
EcoRI	0.5	
XhoI	0.5	
	50 μ l.	

GET Paperwork ready for

HTPAN08504 + PQE6
 HSKBN09 + PQE7

HMSAF22 + GPPA2
 HMSAF22 3' Δ + GPPA2

To Submit to Protein Expression.

Make Primers to Clone into pCDNA
3' HA Tag

HE2PM21 - Soluble VEGF

HMSAF22

HT4SB02 } Possible Secreted Protein

HSKBN09

Make Primers to Clone into PQE6

HTPAN08504 - at New ATG Start

1st for Cloning Cre facility

Inoculate for max pups

HTPAN08 + PQE6 } TB + Amp Kan

HSKBN09 + PQE7

HMSAF22 + GPPA2

HMSAF22 3' Δ + GPPA2 } TB + Amp

2/7/96

Spin through G-25 Column - 1.3K.
Collect fluids through
Count tube

	SAM	POS	CH	CPM	2SIGZ	TIME	EL TIME	AVG H#	RCHZ
NSong	HSLED86	1	200	1 1258393.25	0.46	0.15	1.58	61.0 62.9×10^6	0.00
NSong	HSBV83	2	201	1 1640159.88	0.40	0.15	3.38	112.0 82×10^6	0.00
NSong	HCUER32.3	202	1	1542533.25	0.42	0.15	5.16	87.0 77.15×10^6	0.00

Specific Activity cpm/ug

HSLED86 1.3×10^9 cpm/ug

HSBV83 1.4×10^9 cpm/ug

HCUER32 1.5×10^9 cpm/ug

Give Probes to Mark Porter on 2nd floor

Diagnose Max Preps

Spin Cultures 5K 15min
Make Glycerol Stock of Cultures
Store -80°C

Pour off Supernatant

Resuspend pellet 10ml Buffer P1

Let sit at room Temp 10min

Add 10ml of Buffer P2

Add 10ml of ice cold Buffer P3

Incubate on ice 30min

Spin 8K 30min at 4°C

Equilibrate a tip 500 with 10ml

Buffer QBT

Apply supernatant to a equilibrated

Column through a Kim wipe

Allow Supernatant to completely

flow through Column Bed.

2/7/96


Wash column 2x with 30 ml
Buffer QC
Elute in 15 ml Buffer QC
Add 0.7 volumes (10.50 ml) of isopropanol
Mix well
Spin 8K 30 min
Pour off supernatant
Wash pellet 15 ml 70% ethanol
Spin 8K 15 min
Pour off supernatant
Allow pellet to air dry at
Room temp O/N

2/8/96

Diagen MaxIS Continued
Resuspend pellets in 200 μ l
TE
Transfer to fresh tube
Read OD 260/280

Sample ID	abs		260.0 nm		280.0 nm	
	260.0 nm	280.0 nm	260.0 nm	280.0 nm	260.0 nm	280.0 nm
1) HMSAF22	0.0077	0.0040	2.4654	0.4055	0.1 μ g/ μ l	
2) HMSAF22	0.0563	0.0360	1.5652	0.6389	0.56 μ g/ μ l	
3) HSEBAC9	0.2331	0.1531	1.5226	0.6568	2.33 μ g/ μ l	
4) HTPAN08	0.0790	0.0507	1.5575	0.6421	0.79 μ g/ μ l	

Run gel on gel with 1 kb ladder

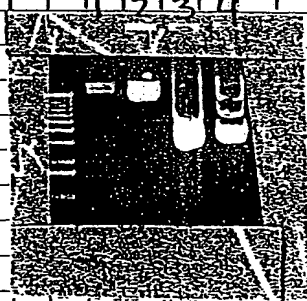

1) HMSAF22 + GPPA2
2) HMSAF22 3' + GPPA2
3) HSEBAC9 + PDE2
4) HTPAN08 + PDE6

2/8/96

Set-up Digestions


	DNA	10X	H ₂ O	Enzs
1 HMSAF22 + GPPAZ	10	2	7.6	0.2 Bgl II + Asp7
2 HMSAF22 3Δ + GPPAZ	1.8	2	15.8	0.2 Bgl II + Asp7
3 HSKONOR + PQE7	0.4	2	17.2	0.2 Sap I + Hinf
4 HTPANOSH + PQE7	1.3	2	16.3	0.2 Xba I / Hinf

Incubate 37°C 4 hrs

Run 10 μl on gel w/ 1 kb ladder

#2 + #4 look good + correct

#3? more than one sight?

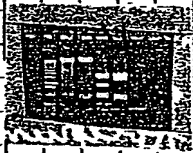Incubate 37°C O/W.
Run again tomorrowSet up reactions for Pfu using
PCR Optimization Buffers
from Stratagene

10 mM Tris-HCl	MgCl ₂	25 mM KCl	75 mM KCl
pH 8.3	1.5 mM	Opti-Prime™ 1x buffer #1	Opti-Prime™ 1x buffer #2
pH 8.3	3.5 mM	Opti-Prime™ 1x buffer #3	Opti-Prime™ 1x buffer #4
pH 8.8	1.5 mM	Opti-Prime™ 1x buffer #5	Opti-Prime™ 1x buffer #6
pH 8.8	3.5 mM	Opti-Prime™ 1x buffer #7	Opti-Prime™ 1x buffer #8
pH 9.2	1.5 mM	Opti-Prime™ 1x buffer #9	Opti-Prime™ 1x buffer #10
pH 9.2	3.5 mM	Opti-Prime™ 1x buffer #11	Opti-Prime™ 1x buffer #12

10

2/8/96

Use Different Sized DNA inserts

1	HIPB411515	5' Bam	3' Xba
2	HIPAN07504	5' Nco	3' Xho
3	HE8CT26	5' Bam	3' Asp
4	HT3AB35	5' Bam	3' Hinf

Use primers at - 20 pmol / reaction
 or 2 µg / reaction
 DNA - at 100 ng / µl

Set up 3 Reaction Tubes per DNA sample

Add 5 µl of 10x Buffers 1-12

add 5 µl of Pfu Buffer to #13

- Make Cocktail in the following Order

* Make sure all reagents + Tubes are on ice *

	1	2	3	4
H ₂ O	41.2	270.6	231.2	469.6
50x Master Mix	12.5	12.5	12.5	12.5
10x dNTP	60	60	60	60
5' Primer	67.4	67.4	92.6	10
3' Primer	5.6	147	161.2	5.4
DNA	1	1	1	1
Pfu	6.5	6.5	6.5	6.5
	56.5	56.5	56.5	56.5

45 µl / Reaction Tubes


Aliquot into tubes just before ready
 to start PCR

2/8/96

PCR PROG #1.

95°C 5min
 95°C 1min
 65.5°C 1min
 72°C 1min + 5 sec.
 72°C 7min
 4°C Hold.

Run 10 μ l on gel w/ 11 kb ladder.

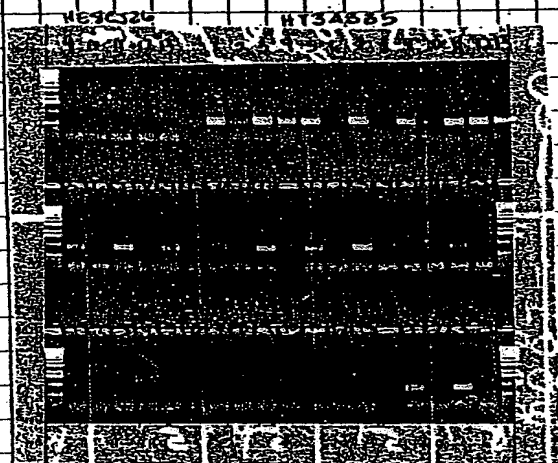

looks like HTPB411S15
 was too large
 to amplify -
 Need longer
 extension time

For the most part
 it looks like the
 longer fragments
 like odd # Buffers
 #1, 3, 5, 7.

Try HTPB411S15 Again
 using a longer
 extension time?

2/9/96

Set up Reactions for Core
 Cloning

Total of 8 samples -

6 mine

2 from John Green

2/9/96

			5' Primer	3' Primer
A	HUSAQ05	10 ng/ul	4624	4623
B	HNEDW90	10 ng/ul	4626	4625
C	HE8C326	250 ng/ul	4602	14603
D	H.T4SB02	250 ng/ul	4636	4637
E	H.MSAF22	250 ng/ul	4634	4638
F	H.E2PM211	250 ng/ul	4633	4629
G	H.SKBN09	250 ng/ul	4635	4630
H	H.TPAND08	250 ng/ul	4632	14388

Set-up reactions using Buffers
1, 3, 5, 7, 9
Total of 200 ul of each Reaction

OK ICE *

	A	B	C	D	E	F	G	H
H ₂ O	75.3	75.2	67.2	77.9	77.8	77.7	77.7	75.7
70x M.M.	20	20	20	20	20	20	20	20
10x dNTP	100	100	100	100	100	100	100	100
5' Primer	6	7	37	17	12	14	19	13
3' Primer	6	6	65	14	2	9.9	1.1	20.0
DNA	10	10	1	1	1	1	1	1
PFU (2.5%)	5	5	5	5	5	5	5	5
	900 ul	900 ul	900 ul	900 ul	900 ul	900 ul	900 ul	900 ul

Aliquot 90 ul into each tube

PCR:




95°C 5 min
95°C 1 min
55°C 1 min
72°C 2 min
72°C 7.5 min
4°C Hold

5 4629 0.3507 0.1821
6 4630 0.2808 0.1819
7 4631 0.2172 0.1518
8 4632 0.2431 0.1517
9 4633 0.2173 0.1450
10 4634 0.2653 0.1770
11 4635 0.1585 0.0989
12 4636 0.1778 0.1265

1.8258 0.5477 2.51
1.5434 0.6479 1.85
1.4310 0.6988 1.43
1.6032 0.6238 1.46
1.4989 0.6672 1.43
1.4992 0.6670 1.75
1.6128 0.6201 1.05
1.4057 0.7114 1.17

2/9/96

Run 5ul on gel w/ 1 Kbladder

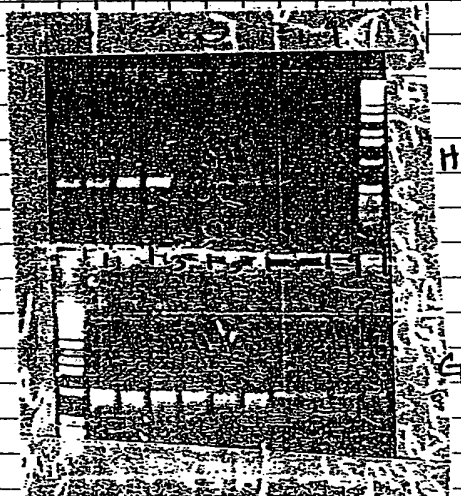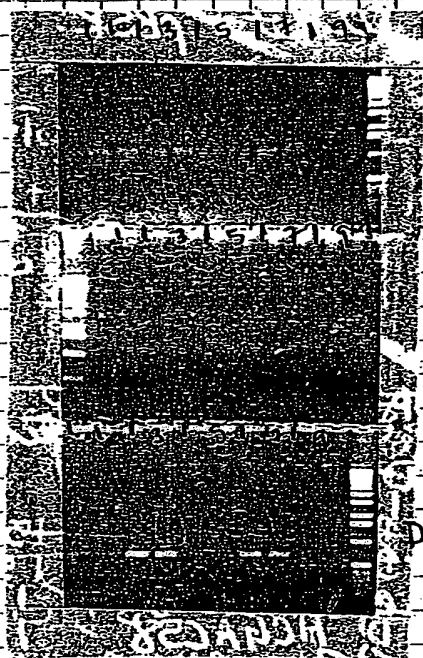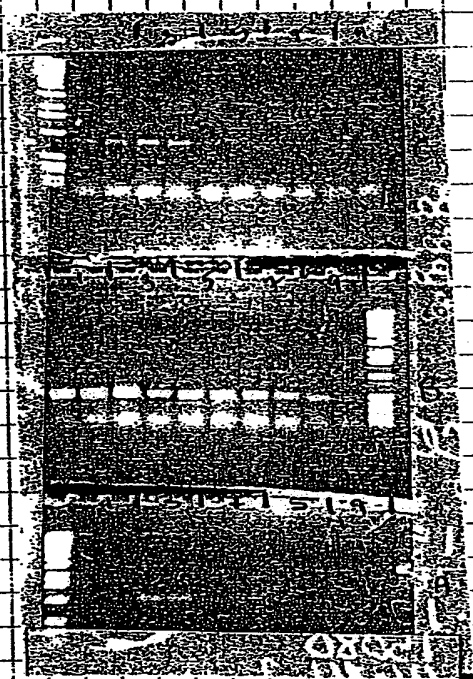

Combine Rixms that
Worked into 1-2 tubes.

Ppt w/ equal Volumes
13% PEG ~~DETA~~ 1.4M NaCl
Vortex well to mix
Spin 15 min
Pour off Supernatant

Store -20°C till
Monday

14

2/9/96

ROXANNE DUAN

Brought up samples for Cloning

Digest pBSK

DNA	3.7
10x Bam	5
H ₂ O	40.3
BamHI	1
	50

Incubate 37°C


2/12/96

Set-up Reactions for Cloning

0.8kb	A)	HSGSA61	15091	15077
0.8kb	B)	HSGSA61	15091	15090
0.6kb	C)	HCGAC58	15092	15089
0.6kb	D)	HCGAC58	15092	15080
0.5kb	E)	HAVCC34	4639	4641

	A	B	C	D	E
di H ₂ O	734.8	734.8	734.8	734.8	734.8
50XMM	20	20	20	20	20
10x dNTP	100	100	100	100	100
5' Primer	20	20	20	20	20
3' Primer	20	20	20	20	20
DNA	0.2	0.2	0.2	0.2	0.2
PEU	5	5	5	5	5
	900	900	900	900	900

2/12/96

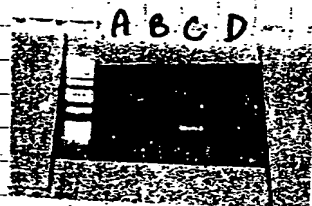

Set-up Digests.

A	HSGSA61	5091 + 15077	Bam → PBL
B	HSGSA61	5091 + 15090	Bam → PBL
C	HCGAC58	15092 + 15089	Bam → PBL
D	HCGAC58	15092 + 15080	Bam/Xho → pcDNA-3'HA

	A	B	C	D
DNA	20	20	20	20
10X	5	5	5	5
H ₂ O	24	24	24	23
Enz I	1	1	1	1
Enz 2	—	—	—	1
Incubate	37°C	37°C	37°C	50°C

Spin Samples from 2/9/96.

15 min
 Pour off Supernatant
 Wash pellet 1000 ul 70% ethanol
 Spin 5 min
 Pour off Supernatant - let air dry 15 min
 Resuspend pellet in 50 ul TE

300 bp Digest

Run gel on gel w/ 1 kb ladder



1	HUS4005	1.12
2	HWFDW90	0.48
3	HEGCT26	0.7
4	HT4SBD2	0.65
5	HMSAE22	0.92
6	HE2PM21	1.1 Kb
7	HSLBN09	1 Kb
8	HTPON26	0.57

500 ul
 Digest
 pgs 20

(pg 16)

2/12/96

A	HUSAQ05	Bgl II / xho I	#2
B	HNFDW10	Bam HI / xho I	Bam
C	HE8J26	Bam I / Asp.	B
D	HT4S602	Bam / xho	Bam
*E	HMSAF22	Bcl / xho	M
F	HEZPM21	Bam / xho	Bam
*G	HSKBND9	Bcl / xho	M
*H	HTPAN08	Bsp HI / HLL	B

	A	B	C	D	E	F	G	H
DNA	20	10	20	20	30	20	20	20
10X	5	5	5	5	5	5	5	5
H ₂ O	23	33	23	23	13	23	23	23
Enz 1	1	1	1	1	1	1	1	1
Enz 2	1	1	1	1	1	1	1	1

- * For Bcl I Digests - Digest first in xho I at 37°C - then Add Bcl and incubate 50°C
- For Bsp HI - Need to use isosomer Rea. I - No Bsp HI around

incubate at 37°C. 6 hrs

Aliquoted HIPAN08 SOL 4/3/12



for Patent

25 Tubes w/ 100ng DNT at 10ng/ul

Give to Steve Ruben

2/13/96

Run Digests on LMP Gels
 80V 1 1/2 hrs.
 Cut out Gel Slice
 Take Pictures

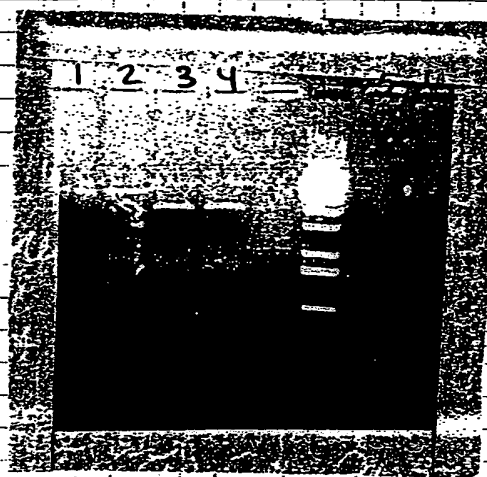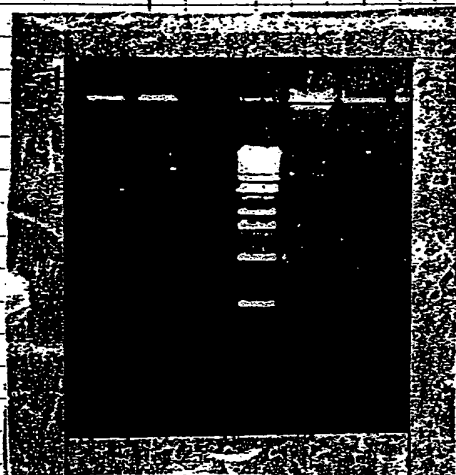

1 HCGAC58	Bam/Xho	pCDNA 3' HA
2 HCGAC58	Bam	pD10
3 HSGSA61	Bam	pD10
4 HSGSA61	Bam	pBSK
5 HNF DW90	Bam/Xho	pCDNA 3' HA
6 HUSAQ05	Bcl I/Xho	pCDNA 3' HA
7 HEZPM21	Bam/Kpn	pCDNA 3' HA
8 HE8BT26	Bam/Ksp	pA2
9 HMSAFC2	Bcl I/Xho	pCDNA 3' HA
10 HSCBNO9	Bcl I/Xho	pCDNA 3' HA
11 HTY5602	Bam/Xho	pCDNA 3' HA
12 HTPAN08	BspH I/HuI	pQEG

Gene Clean fragment

Set up ligations

2/13/90

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
1 HCGAC58 Bam/Kho	10															
2 HCGAC58 Bam		10														
3 HSGSA161 Bam			10													
5 HNFDA90 Bam/Kho				10												
6 HUSAQ05 Bsp/Kho					10											
7 HEZFM21 Bam/Kho						10										
8 HEXCJ26 Bam/Asp							10									
9 HMSAF22 Bsp/Kho								10								
10 HSKBN09 Bcl/Kho									10							
11 HTUS802 Bam/Kho										10						
12 HTPAN08 Bsp/HII											10					
PCDNA 3' HA Bam/Kho	0.5			0.5	0.5	0.5		0.5	0.5	0.5		0.5				
PD10 Bam		0.5	0.5										0.5			
PA2 Bam/ASP							0.5							0.5		
PQE6 Nco/HII											0.5				0.5	
10X T4 Buffer	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2
T4 Ligase	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
H ₂ O	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
Total	20	20	20	20	20	20	20	20	20	20	20	20	20	20	20	20

incubate 16°C overnight

2/13/96

Inoculate LBT Amp with:


HIPANOS04

GPPA2⁴

pCDNA 3'HA

Incubate overnight 37°C w/aeration
overnight

2/14/96

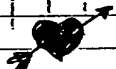

VALENTINE'S DAY

Qiagen Maxi Prep

HIPANOS04

GPPA2

PCDNA 3'HA

Spin 5K 15 min

Pour off Supernatant

Resuspend pellet 10 ml P1 Buffer

Add 10 ml P2 Buffer

Add 10 ml P3 Buffer

Incubate on ice 30 min

Spin 8K 30 min

Transfer Supernatant through
Ken tube to equilibrated
tip 500

Allow to flow through

Wash tip 3x with 30 ml

P1 Buffer

Elute with 15 ml Buffer QF

Collect - add 0.7 volumes (10.5 ml)

of isopropanol

Spin 8K 30 min

Pour off Supernatant

Wash pellet 15 ml 70% Ethanol

Spin 8K 15 min

25

2/14/96

Pour off Supernatant
 Allow pellet to air dry
 Resuspend pellet in 200 μ l TE
 Read OD₂₆₀/280

Sample ID	abs 260.0 nm	abs 280.0 nm	260.0 nm 280.0 nm	280.0 nm 260.0 nm	
1	0.0453	0.0388	1.2477	0.8015	0.46 μ g/ μ l
2 HTPANOSY	0.0457	0.0366	1.2477	0.8015	1.54 μ g/ μ l
3 GIPPAZ	0.1593	0.1044	1.5253	0.6556	0.46 μ g/ μ l
4 pCDNA 3' HA	0.0458	0.0320	1.4240	0.7022	

Set-up Digestions of the pC4. Prom
 2/13.

PC4 - 0.96 μ g/ μ l.

pCDNA 3'

	BamHI	XbaI	Asp718	BamHI/XbaI	Bam/Xba
DNA	5.21	5.21	5.21	5.21	10.9
10X Buffer	5 Bam	5 #2	5 8	5 Bam	5
H ₂ O	37.79	37.79	37.79	36.77	30.1
Enzyme 1	1	1	1	1	1
Enzyme 2	1	1	1	1	1
	50	50	50	50	50

incubate 37°C O/N

Transform ligations into Chemically
 Competent Cells.
 for pCDNA + PAZ - DH5 α
 pOE + PD10 - M15 up 5

2/14/96

Thaw Cells on ice
 aliquot 100ul into prechilled
 tubes
 add 10ul of ligation mix
 incubate on ice 7 hours
 heat 42°C 45 sec
 chill on ice
 add 400ul LB
 incubate 37°C 1hr
 plate into LB plates
 pCDNA + PAZ - LB + Amp
 pAE + PD10 - LB + Amp / Kan
 incubate 37°C overnight
 use (+) control PAZ idig
 & pAE 100 long
 use (-) control of cells only.

2/15/96

Pick Colonies into LB + Amp
 or LB + Amp / Kan in 96 well Dish
 incubate 37°C w/ aeration 4 hrs.
 Set up PCR

	A	B	C	D	E	F	G	H	I	J	K
10xPCR	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2
10xDNTP	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2	3.2
H ₂ O	21.4	22.3	22.3	20.4	21.4	20.3	21.4	19.4	20.4	21.4	28.2
Tag	0.2	10.2	0.2	0.2	0.2	0.2	0.2	0.2	0.2	0.2	0.2
Cell	2	2	2	2	2	2	2	2	2	2	2
Primer #1	T71	PD10	PD10	T71	T71	PAZa1	T71	T71	T71	T71	PD10
Primer #2	1	1	1	2	1	3	1	3	2	1	01


2/15/96

Primers

A	HCGAG58	pCDNA	#15000	10pmol/ μ l	+ T7
B	HCGAG58	PD10	#15000	(8pmol) μ l	+ PGE9/10
C	HSGSAG1	PD10	#14810	10pmol/ μ l	+ PGE9/10
D	HNF DW90	pCDNA	#11693	7.7pmol/ μ l	+ T7
E	HUSAQ05	pCDNA	#14118	14.4pmol/ μ l	+ T7
F	HE2PM21	pCDNA	FPO2	3.2pmol/ μ l	+ T7
G	HESCS26	PA2	FPO4 #1177	(10.5pmol) μ l	+ T7 & 20
H	HMSAF22	pCDNA	FPO7	11577 (3.2pmol/ μ l)	+ T7
I	HTSKB09	pCDNA	FPO6	11463 (5.6pmol/ μ l)	+ T7
J	HTYS02	pCDNA	FPO5	11746 (7.0pmol) μ l	+ T7
K	HTAN08	PDEG	FPI4		+ PGE9/10

PCR Program lab.

Run 10 min on qpt w / 1 to ladder

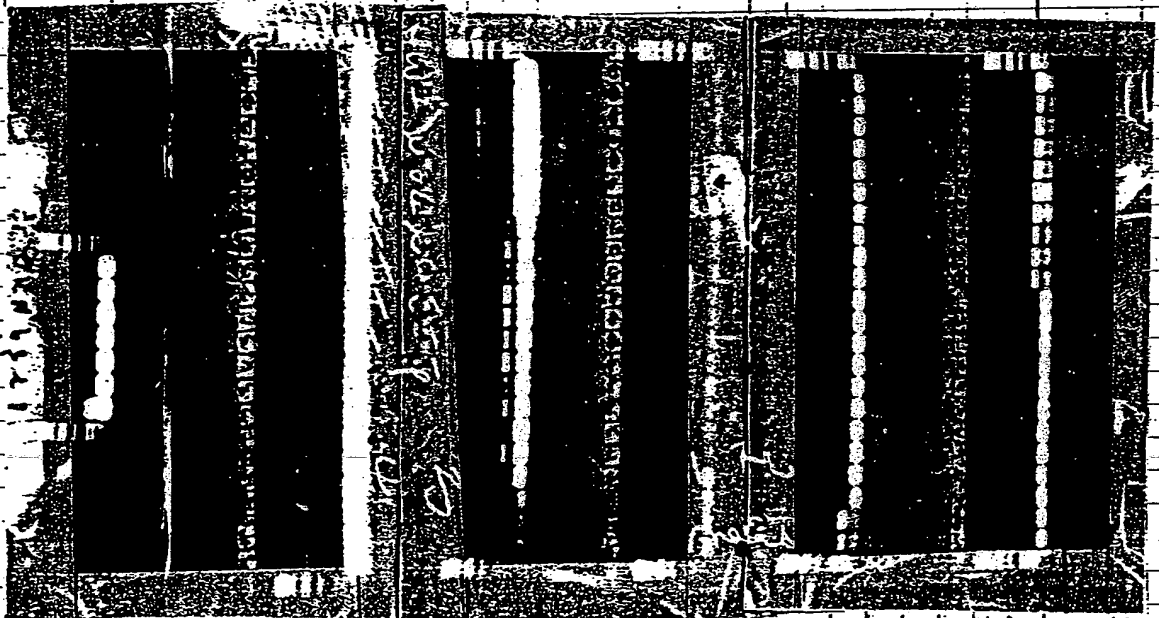

2/15/96

Inoculate for mini Prips

Try to PCR insert with Tag

Take some of PCR Fragment Generated from 2/9 & 2/12/96.

PCR using same primers.
+ 0.5 μ l Tag. - Total of 500 μ l Reactions.

PCR.

95°C	1/3 5min	} 25x
95°C	1min	
55°C	1min	
72°C	1.45min	
72°C	7:30min	

Run 5ul on gel w/ 1Kb ladder



- 1 HCGAC58
- 2 HCGAC58
- 3 HSGSAG1
- 4 HSGSAG1
- 5 HNFQ390
- 6 HUSAQ05
- 7 HE2PM21
- 8 HEC726
- 9 HMSAF22
- 10 HSKR109
- 11 HTUS002
- 12 HTPAN08

PEG PPT.

Spin pellet 1000^g for 5min of 30% ethanol.
Pour off supernatant

2/15/96


Allow pellet to Air Dry
 Resuspend in 100 μ l TE
 Set up Digests

	1	2	3	4	5	6	7	8	9	10	11	12
DNA	25	25	25	25	25	25	25	25	25	25	25	25
10X Buffer	5	5	5	5	5	5	5	5	5	5	5	5
H ₂ O	18	19	19	19	18	18	18	18	18	18	18	18
Enz 1	1 Xho	1 Bam	1 Bam	1 Bam	1 Bam	1 BglII	1 Bam	1 Bam	1 Bcl	1 Bcl	1 Bam	1 Bcl
Enz 2	1 Bam				1 Xho	1 Xho	1 Asp	1 Xho	1 Xho	1 Xho	1 Xho	1 Hcl
	50	50	50	50	50	50	50	50	50	50	50	50

Incubate 37°C O/N.

W/ Bcl Digests - incubate at 50°C for 4 hrs before Running on Low P Gel

Run 2 μ l of PC4 Digests mg/w/1 kb ladder

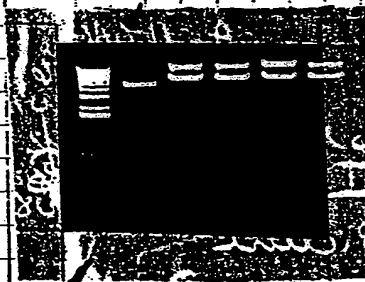


1 pCDNA 3'HA Bam/XhoI
 2 Asp 718
 3 Bam #1
 4 Xba
 5 Bam/Xba } pC4

Run 1/2 of #1 & 5
 on 1% P Gel
 with PC1 Bam/Xba
 and 1 kb ladder

Run at 80V 2 hrs
 Cut out DNA Fragment
 Take P. Clunk

30

2/15/96

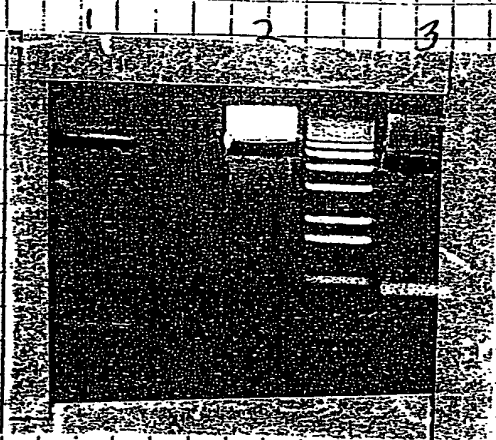

1 PCL Bam/Xba
 2 pCH Bam/Xba
 3 pCDNA 3 HA
 Bam/Xba


Store -20°C
 O/N.

Sample ID	abs 260.0 nm	abs 280.0 nm		260.0 nm	280.0 nm
				280.0 nm	260.0 nm
1 PCL 5' (4)	0.1281	0.0779	24 mer	1.8454	0.6078
2 MCS rppi 3' (10)	0.1530	0.1037	21 mer	1.4745	0.6782

0.85 µg/µl 106.8 pmol/µl
 1 µg/µl 145.7 pmol/µl

2/16/96

Run Digests on 0.8% LMP Gel w/
 1 kb ladder
 Run 80V 2hrs
 Cut out gel slices & take picture

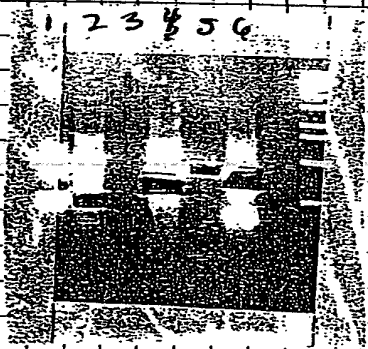

2/16/96

Melt Gel slices at 68°C
Set up ligations

	A	B	C	D	E	F	G	H	I	J	K	L
2 HCLAKES B	10											
4 HSCGSA61 B		10										
5 HNFOW 90 O/X			10									
6 HNSAQOS B/X				10								
7 HEP2M21 B/X					10							
8 HSCG26 B/X						10						
9 HNSAFU B/X							10					
10 HSEBNO9 B/X								10				
11 HTY5B04 B/X									10			
12 HTPAN28 B/HW										10		
10X Liga Buf.	5	5	5	5	5	5	5	5	5	5	5	5
H ₂ O	23.5	23.5	23.5	23.5	23.5	23.5	23.5	23.5	23.5	23.5	23.5	23.5
T4 Ligase	1	1	1	1	1	1	1	1	1	1	1	1
pcDNA 3.1 HAB4			0.5	0.5	0.5		0.5	0.5	0.5			
pDIO Barn	0.5											
pBSK Barn		0.5										
pAZ Barn / ksp						0.5						
pQEG N/HW										0.5		
	50	50	50	50	50	50	50	90	50	50	50	50

Incubate at Room Temp 2 hrs.
Store 4°C
Company Closed 1pm due to weather.

Baling Manipulations


Spin 2ml Culture
Remove Supernatant
Resuspend pellet in 700ul STET +
Wage
Heat to 100°C 5min

2/16/96

Spin 15min
 Transfer 650 μ l of Supernatant to
 Add Fresh Tube
 Vortex 650 μ l of 13% PEG / 1.6M NaCl


Spin 15min
 Pour off Supernatant
 Wash Pellet 1000 μ l 70% Ethanol
 Spin 5min
 Remove Supernatant
 Allow Pellet to Air Dry 15min
 Resuspend in 100 μ l THF

Run 1 μ l on gel w/ 1 kb ladder

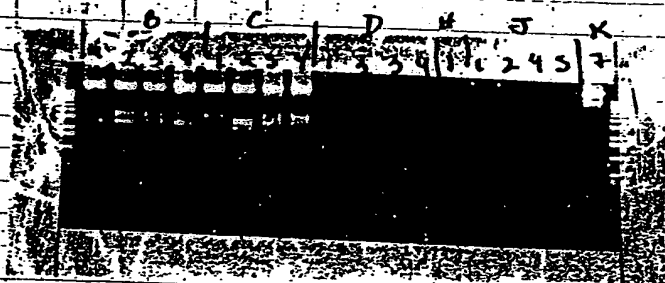

B HCGAG58-PD10
 C H8SA61-PD10
 D HNFDF10-PC3HA
 J H4SAF22-PC3HA
 K HTPAN88-PD10

Set up digestion B/C - with Bam
 and K w/ Hind III. ~~at 37°C~~ 37°C

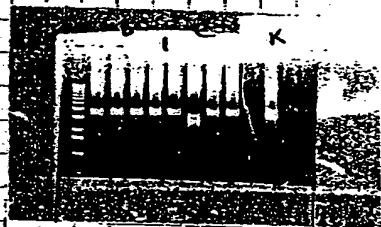

None Digested

Redo Min 5 on
 Monday

2/19/96 OFF President's Day

2/19/96

Transform ligations from 2/16/96

Transform PQE & PD10 constructs into
M15 up4

Transform PAZ, PBSK + pCDNA 3' HA into
DH5 α cells.

Thaw Chemically Competent cells
on ice

Aliquot 100 μ l into Sterile Tubes

Add 10 μ l of Ligation Reaction

Incubate on ice 1 hr

Heat 42 $^{\circ}$ C for 45 Sec

Quick chill on ice

Add 400 μ l LB

Incubate 37 $^{\circ}$ C 1 hr

plate onto LB + Antibiotics plates.

M15 up4 cells - LB + Amp/Kan

DH5 α cells - LB + Amp

Incubate 37 $^{\circ}$ C O/N

Inoculate for mini preps - 2/16/96

2/20/96

plates are again contaminated

- just LB + Amp

picked clones anyway

into 96 well Dish


incubate 37 $^{\circ}$ C O/N

mini preps -

Do Promega preps

2/20/96

Spin 2 ml culture 2 min
 Remove Supernatant
 Resuspend pellet in 500 μ l Resuspension Buffer
 Add 500 μ l Cell Lysis Buffer
 Mix by inverting tube several times
 Add 500 μ l of Neutralization Solution - mix well
 Spin 15 min
 Transfer Supernatant to fresh tube with 500 μ l of Resin Slurry
 incubate 5 min at RT
 Use Cologun with 3CC syringe
 Attach
 Apply Resin & Supernatant to Column
 Vacuum through
 Wash 2 x 1 ml Wash Buffer
 Spin Column 1 min
 Apply 150 μ l TE - heated to 68°C
 Let sit 5 min
 Spin 1 min & collect
 Run Gel on gel w/ 1 kb ladder

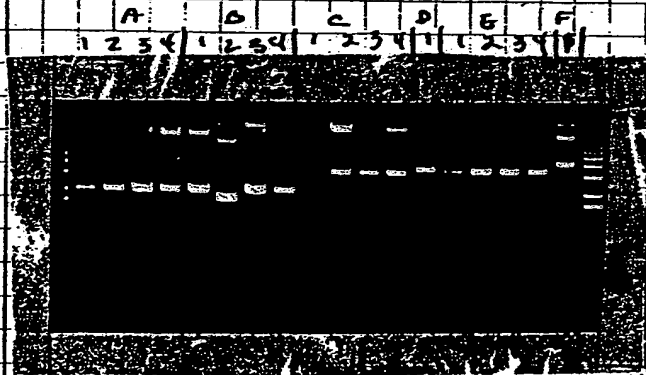

A - HCGAG58 + PD10
 B - HSGSAG1 + PD10
 C - HNF200 + PC3
 D - HUSAF22 + PC3
 E - HT45802 + PC3
 F - HTPAN08 + PC3

PCR Cultures

2/21/96

A - HCGACS8 Bam + PD10 Bam
 C - HNF1W90 Bam/Xho + pCDNA 3'HA
 D - HUSAQ05 BspEI/Xho + pCDNA 3'HA
 E - HE2pm21 Bam/Xho + pCDNA 3'HA
 F - HE8CJ26 Bam/Asp + PAZ
 G - Hm5AF22 Bcl/Xho + pCDNA 3'HA
 H - HSKBND9 Bcl/Xho + pCDNA 3'HA
 I - HT4SB02 Bam/Xho + pCDNA 3'HA
 J - HTPAN08 BspH/HII + pDEG0

C, D, E, G, I, J. - pCDNA 3'HA

T7	1	72x
Sp6	2	72
10x dNTP	3.2	144
10x PCR	3.2	230.4
H ₂ O	20.4	230.4
Cult	2	1468.8
Tag	0.2	14.4
	30ul	30ul/Tube

E - PAZ		12x
T7 Bsp	0.05	0.6
3 Asp 718	3	36
10x dNTP	3.2	38.4
10x PCR	3.2	38.4
H ₂ O	20.4	244.8
Cult	2	
Tag	0.2	2.4
	32	30.2/tube

36


2/21/96

A.5 - PD10 + POE60.

3 POE	0.2	3.3X
17 POE	1	6.6
10X	3.2	33
10X	3.2	105.6
H2O	20.2	105.6
Cult	2	732.6
Lat	0.2	6.6
	30ul	30ul/well

PCR Program (do)

Run 10ul on gel w/ 11kb ladder

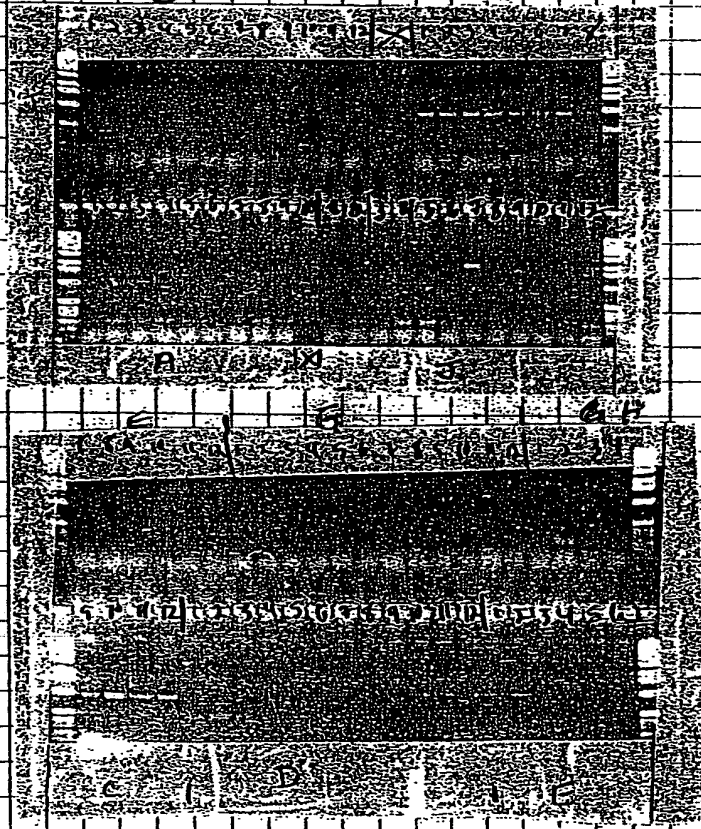

2/21/96



- Unincubate some pBSK - HSGSACI Bam.
into Sm. TB + Amp.

- Unincubate
HNE DW90 CI-4
HUSAQ05 DI-DS.
HEZPM2 E2.

into TB + Amp. for E. coli Mini prep


2/22/96

Proteoglycan Mini prep DNA.
See pg 33-34 for Protocol.

2/24/96

Setup

Run 1ul on gel with 1 kb ladder

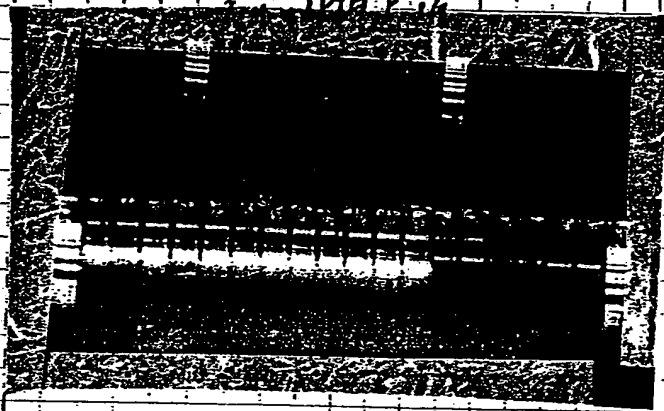

A - HSGSAGI + PBSK
 B - HNF-1α
 C - HUSARα5
 D - HEZPM21
 E - H' MSAF22
 F - HSEBNO9
 G - H7PAND8

Sett - ca Digests

pCDNA 3.7HA			pBSK		
HIII / Xba			Xba I		
		(22X)			(13X)
DNA	10		DNA	10	
100X IX#2	3	66	100X IX#2	3	39
BSA	0.3	6.6	BSA (100X)	0.3	39
H ₂ O	16.3	38.6	H ₂ O	22.5	49.5
HIII	0.2	4.4	Xba I	0.2	8.6
Xba I	0.2	4.4		30	36ul / tube
	30	32ul / tube			

pDIO + pDGO

Eco RI / HIII

		12.4	Incubate 37°C 4hr	
DNA	10			
100X IX#2	3	36		
BSA	0.3	3.6		
H ₂ O	16.3	195.6		
HIII	0.2	3.6		
Eco RI	0.2	3.6		
	30ul	30ul / tube		

2/27/96

Set up Digest.

- Eco RI / Xho I.

Digest DNA to 250 ng / μ l

DNA	8	5x
10x EcoRI	5	25
H ₂ O	36	180
Eco RI	0.5	2.5
Xho I	0.5	2.5
	50	40 μ l / Tube

also Digest HCE 36226 1.32 μ g / μ l
HCE 5F84 1.57 μ g / μ l

DNA	4
10x EcoRI	5
H ₂ O	40
Eco RI	0.5
Xho I	0.5
	50.0

Incubate 37°C O/N.

Cassie + I Packed up Lab #15 on
3rd floor 9:20. will be
moving to Lab #5 2nd floor
9:20. J

2/28/96

Cassie + I Spent the whole Day
moving Boxes & Moving into
New Space!

2/29/96

Set-up PCR of colonies picked
2/26/96

- | | |
|---------------------------|--------------------------------|
| 1. HMSAE22 + GPPAZ | T7 + 3' Bof II |
| 2. HMSAE22 3'Δ + GPPAZ | 17 Bacc + 3'Δ A _{top} |
| 3. HCEED20 + GPPAZ | 13627 + 907 T7 Bacc |
| 4. HCEG495 + pCDNA 3' HA | } cho 3' + 5' |
| 5. HFGAm58 + pCDNA 3' HA | |
| 6. HATCK 87 + pCDNA 3' HA | |

	0.1	0.6
T7 Bacc	0.1	0.6
3' Bof II	2	12
10x dNTP	3.2	19.2
10x PCR	3.2	19.2
H ₂ O	21.2	127.2
Taq	0.3	1.8
Cit. 1	2	
	32ul	30ul/tube

	0.1	1.2
T7 Bacc	0.1	1.2
3'Δ A _{top}	3	36
10x dNTP	3.2	38.4
10x PCR	3.2	38.4
H ₂ O	20.2	242.4
Taq	0.3	3.6
Cit. 1	2	
	30.0	30ul/tube

	0.1	25x
T7	0.1	2.5
3' 13627	0.4	30
10x dNTP	3.2	80
10x PCR	3.2	80
H ₂ O	21.8	545
Taq	0.3	7.5
Cit. 1	2	
	32ul	30ul/tube

	0.02	75x
3' cho	0.02	1.5
5' cho	0.02	1.5
10x dNTP	3.2	270
10x PCR	3.2	270
H ₂ O	23.26	174.5
Taq	0.3	22.5
Cit. 1	2	
	32ul	30ul/tube

PCR Prog	66.
95°C	5min
95°C	30sec
55°C	30sec
72°C	1min
72°C	10min
4°C	hold

30x

2/28/96

3/5/96

Count rel.

	SAM	PUS	CH	CPM	2S16%	TIME	EL TIME	AVG H#
HCE3026	1	296	1	728784.00	0.47	0.25	1.72	60.0
HME604	2	297	1	815745.00	0.30	0.20	3.48	61.0
HCE588	3	298	1	890950.00	0.47	0.20	5.26	65.0

7.3×10^5 /ul
 8.2×10^5 /ul
 8.7×10^5 /ul

Give Probe + DNA strips to
Main Poster.

Phenol Chloroform Extract Gel slides
 from 2/16/96 = Pg 30
 from 2/29/96 = Pg 62.

Add TE to 500ul

Heat 68°C

Add 500ul heated phenol 10

Mix well

Spin 10min

Transfer Supernatant to fresh tube

Add 500ul phenol

Mix well

Spin 10min

Transfer Supernatant to fresh tube

Add 800ul SFA

Mix well

Spin 10min

Transfer Supernatant to fresh tube

Add 20ul 3M NaAc. pH 5.3

Add 1000ul 100% ethanol


Mix well

Let sit on ice 10min

Spin 15min

3/5/94

Pour off Supernatant
 Wash pellet 100ul 70% Ethanol
 Spin 5 min
 Pour off Supernatant
 Allow pellet to Air Dry 15 min
 Resuspend pellets in 40ul TE
 Run gel mixed w/ 1 kb ladder

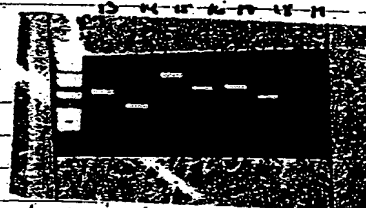

1	HCGAC58	Bam/Vho
2	HCGAC58	Bam
3	HCGAC58 HSGSAC1	Bam/H1
4	HNFDW90	Bam/Vho
5	HUSAQ05	Bgl II/Vho ~800bp
6	HUSAQ05	Bgl II/Vho ~450bp
7	HE2Pma1	Bam/Vho
8	HE8CJ26	Bam/Asp
9	HmSA122	Bcl/Vho
10	HSKBN09	Bcl/Vho
11	HTUS022	Bam/Vho
12	HTPAP08	BspH/HTD
13	HLTB050	EcoRI/Vho
14	HCE3W26	0.7kb Eco/Vho
15	HCE3W26	2.5kb Eco/Vho
16	HTPRD07	Eco/Vho
17	HFTD102	Eco/Vho
18	HPRCA09	Eco/Vho
19	HTXem77	Eco/Vho

3/6/96

Set up ligations of ones that
have not worked for the past

	1	2	3	4	5	6	7	8	9	10	11	12
HCEED20 Asp.	4											
HS45AG1 Bam		4										
HCGAC58 Bam.			4									
HUSAQ BglII/Xho				4								
HUSAQ BglII/Xho					4							
HIPAND8 BspHI/Hinf						4						
HCGAC58 Bam/Sma							4					
PBSK Kpn	0.5							0.5				
PBSK Bam		0.5	0.5						0.5			
PCDNA 3'HA Bam/Xho				0.5	0.5		0.5			0.5		
PDE6 Nco/Hinf						0.5					0.5	
10X Buffer	2	2	2	2	2	2	2	2	2	2	2	2
H ₂ O	12.5	12.5	12.5	12.5	12.5	12.5	12.5	16.5	16.5	16.5	16.5	17
T ₄ Ligase	1	1	1	1	1	1	1	1	1	1	1	1
	20	20	20	20	20	20	20	20	20	20	20	20

ligate @ at RT o/n

Therapeutic Protein Meeting

12, 15 -

Carrice Fischer Presentation
on TNTs

3/7/96

Gathered Plasmid preps of all
Constructs with HG project
code

3/7/96

HTPAN08 - HG03500 - Fas Ligand
 pBL HG03500 - pBSK
 pE9 HG03500 - pD10 System
 AZ HG03500 - PAZ construct
 PHA HG03500 - pCDNA3 HA Tag

HTPB911 - HG06300 - Vellin
 pBL HG06300
 pE9 HG06300
 pE60 HG06300 - pQE60
 AZ HG06300

HCDAAS9 - HG09900 - HE2PM21
 pBL HG09900
 AZ HG09900

HSKBN09 - HG10900
 pBL HG10900
 pE7 HG10900 - pQE7
 AZ HG10900


HT45B02 - HG09700
 pBL HG09700
 AZ HG09700
 PHA HG09700

HMSA002 - HG03700 - HMSAP22
 pBL HG03700

HNBAA26 - HG09800
 pBL HG09800

HNFAA64 - HG10000
 pBL HG10000


3/7/96

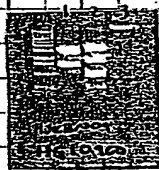

- 1- HSKBN09 Xho/EcoRI
- 2- HSKBN09 + pQE7 - SpeI/HII
- 3- HSKBN09 + PAZ - Bam/Xba



- 1- HT4SB02 Xho/EcoRI
- 2- HT4SB02 + pCDNA3HA Bam/Sma
- 3- HT4SB02 + PAZ Bam/Xba

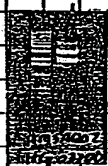

- 1- HMPAF22 Xho/EcoRI



- 1- HMPA79 - Xho/EcoRI
- 2- HMPA79 + pQE60 Bam/BglII

E. coli

Transform *Legatilis* (pg 75)
from 3/6/96.

for pQE6 constructs use M15 sup⁺ cells

for pBSK/pCDNA3HA constructs use XL-1 Blue cells.

Use Chemically Competent frozen 3/7/96
Cells.

- III
- Thaw cells on ice
 - Aliquot 100 μ l into fresh sterile tube
 - Add 10 μ l of ligation mix
 - incubate on ice 1 hour
 - heat 42°C for 45 sec
 - Quick chill on ice
 - Add 100 μ l of LB
 - incubate 37°C 1 hour
 - plate 100 μ l onto plate:
 - PBSK - LB + Amp + IPTG / Xgal
 - pCDNA 3' HA - LB + Amp
 - PQEG - LB + Amp / Kan

Incubate at 37°C O/N.

3/8/96

inoculate 200 μ l of LB + Amp
with: (96 well Dish)
HUSAQ05 + pCDNA 3' HA Tag. (60)
HCGAC58 + PBSK (36)
HCEED20 + PBSK (24)
HSGSAG1 + PBSK (36)

inoculate 200 μ l of LB + Amp / Kan
with: (96 well Dish)
HTPAND8 + PQEG 96.

Incubate all at 37°C w/ aeration

Set up PCRS

03/8/96

HUSAQ05. + PCDNA 3 HA

		62x
T7 PC1	0.3	18.6
14147	1	62
10x dNTP	3.2	198.4
10x PCR	3.2	198.4
H ₂ O	22.2	1386
Taq	0.3	18.6
Cult	2	
	32ul	32ul/tube

HCGAC58 + rcdNA 3 HA

		38x
T7	0.3	11.4
PC 3'	0.1	3.8
10x dNTP	3.2	121.6
10x PCR	3.2	121.6
H ₂ O	22.9	870.2
Taq	0.3	11.4
Cult	2	
	32ul	32ul/tube

HTRAN08 + PQE60.

		100x
T7 PQE	0.02	2
3' PQE	0.25	25
10x dNTP	3.2	320
10x PCR	3.2	320
H ₂ O	23.13	2313
Taq	0.2	20
Cult.	2	
	32ul	32ul/tube

HCEED20 + PBSK

		26x
M13R/E	0.1	2.6
10x dNTP	3.2	83.2
10x PCR	3.2	83.2
H ₂ O	23.3	605.8
Taq	0.2	5.2
Cult.	2	
	32ul	32ul/tube



HSQSA61 + PBSK

		38x
M13R	0.1	13.8
4810 (topred)	1	38
10x dNTP	3.2	121.6
10x PCR	3.2	121.6
H ₂ O	22.3	849.7
Taq	0.2	7.6
Cult.	2	
	32ul	32ul/tube


PCR

95°C	5min
95°C	30sec
55°C	30sec
72°C	1min
72°C	7/2min
4°C	hold



Run 10 of PCR Rxn w/ 1E6 ladder 3/8/96
 HTPANOS + PQEG

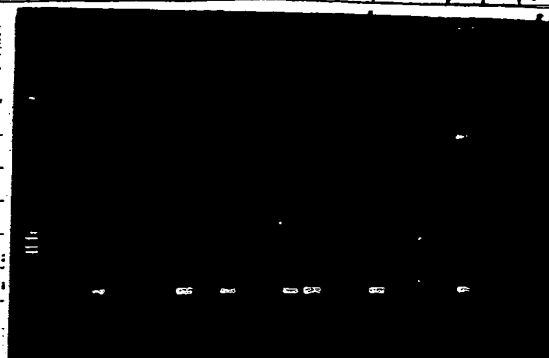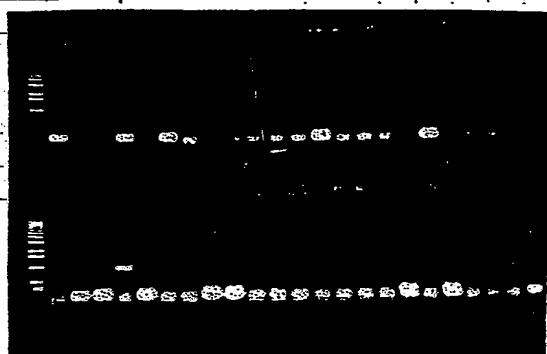

HTPANOS + PQEG

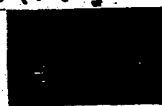

HUSAQOS + PCONA 3' HA

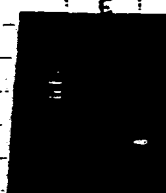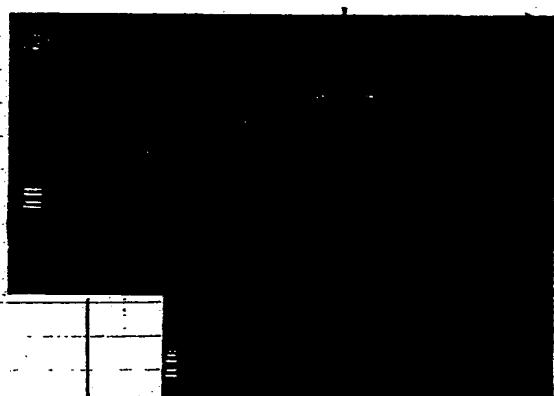

3/15/96

Inoculate TB5 Amp
for mini preps on Monday
leave at RT till Monday

Diagn mini prep

HE2RM21 + 3 HA
HMSAF22 + 3 HA
HSKBN09 + 3 HA
HTUSB02 + 3 HA
HTPAN08 + PDEG

PAT - + Allow to bind till
Monday

Set up ligations

	1	2	3	4	5	6	7
HTAF20 BINCO	5	—	—	—	—	—	—
HTHC008 BamHsp	—	5	—	—	—	—	—
HCEB020 Asp	—	—	5	—	—	—	—
PDEG0 Neo (Bam)	1	—	—	1	—	—	—
PC4 BamHsp	—	1	—	—	1	—	—
PBS1C Asp	—	—	1	—	—	1	—
10X Buffer	2	2	2	2	2	2	2
T4 Ligase	1	1	1	1	1	1	1
H2O	11	11	11	16	16	16	17


Inoculate at RT over weekend

Split TC cells

3/18/96

Diagen Max Cont.

Resuspend PGE construct 400ul TE

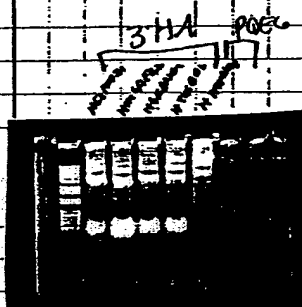Resuspend 3'HA constructs in 200ul TE
Run gel w/ 1kb ladder1 ul of RNA in 3'HA
constructs

qs to 500ul -

Add 500ul PEG/NaCl
mix well

Spin 15min

Wash pellet 1000ul

70% EtOH

Spin 5min

Remove supernatant

Allow pellet to air dry 15min

Resuspend in 100ul TE

Read OD 260 - 280 1.200 Dilution

Sample ID	abs 260.0 nm	abs 280.0 nm	bkg abs 320.0 nm	260.0 nm 280.0 nm	260.0 nm 280.0 nm
1 HE2PM21 + 3'HA	0.1208	0.0775	0.0005	1.5594	0.6413
2 HMEAF22 + 3'HA	0.0934	0.0806	0.0036	1.5765	0.6343
3 HSKBN09 + 3'HA	0.0708	0.0457	0.0044	1.6017	0.6243
4 HT45B02 + 3'HA	0.1050	0.0701	0.0063	1.5474	0.6463
5 HT45B05 + PGE	0.1173	0.0771	0.0073	1.5751	0.6349

Cloned Cell Culture
Set up DigestionHT45B02 2 Com/Max
HE2PM21 3HSKBN09 2 Kpn/Bbs
HMSAF22 3

HT45B02 + PGE 3 Eco/HindIII


3/18/96

DNA	2
10x	3
BSA	0.3
Bam	0.2
Xho	0.2
H ₂ O	24.3
	30 μ l

DNA	2
10x	3
BSA	0.3
Kpn	0.2
Xho	0.2
H ₂ O	24.3
	30 μ l

DNA	2
10x	3
EcoRI	0.2
HindIII	0.2
H ₂ O	24.6
	30 μ l

Incubate 37°C 4hrs
 Run 10 μ l on gel with 1 kb ladder

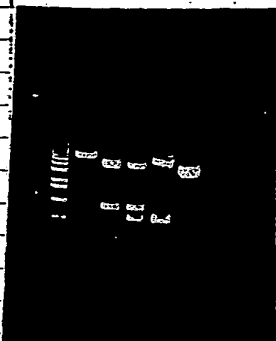

HE2PM2.1 is linearizing
 HMSA F22 - looks correct w/
 RNA contamination
 HSKBN09 - internal seq?
 HT45B02 - incomplete
 digestion

HTPAN08 - incomplete
 digestion

Digest for longer +
 Re Run.

Inoculate TB + Amp or Amp Kan
 with cultures

HSYSA61 Bam + PBSK = A1, A2, B3, B4, D5, D6

HCGAC58 Bam + PBSK = E10, E9, F9, F11, F12
 G1, G2, H3, H4

HTPAN08 (Nco/Hind) + PQE6 = A4, C4, D2

HELBS34 + PCY = D10

HELBS34 + PBSK = G3, G6, G8

HCEED20 + GPPA2 = A1, A2, B3, B4, C5, C6, D7, D8

Incubate 37°C O/N

3/18/96

Transform Ligations into
Chemically Competent cells
XL-1 cells

Chaw cells on ice
Aliquot 100 μ l of cells into fresh tubes
Add 10 μ l of ligations
Incubate ~~on ice~~ on ice 1 hr
Heat 42°C 1 min
Chill on ice
Add 2400 μ l LB
Incubate 37°C 1 hr
Plate onto 2 LB Amp. plates for
PCII + LB Amp. Beta gal for
pBSK.
Incubate 37°C O/N.
Freeze Cole Cells.


3/19/96

Alkaline Lipis Mem Prep

Spin 2ml culture 2 min.
Remove Supernatant
Resuspend pellet 300 μ l Resuspension
Buffer
Add 600 μ l Cell Lysis Buffer
Add 300 μ l Neutralization Buffer
Mix well
Let sit on ice 15 min
Spin 15 min
Transfer Supernatant to fresh tube
Add equal Volume 10% PEG
1.6M NaCl.
Mix well
Spin 15 min
Remove Supernatant

3/19/94

wash pellet 100ul 70% Etanol.
Spin 5min
Remove supernatant
allow pellet to air dry 15min
Resuspend pellets in 150ul TE
Run 2ul on gel with
1k ladder

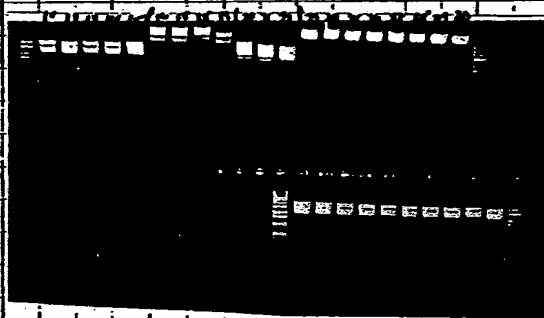1-6 HSGSAG1 +
PBSK7-15 HCGAC58 +
PBSK16-18 HTPAN08 +
PDE619-22 HREHEL05 + Ry
+ PBSK

23-30 HCEED20 + GPPAR

Mini PCR's look good.
Set up Digestions

HSGSAG1 + PBSK
HCGAC58 + PBSK
PBSK

HTPAN08 + PDE6
PDE6

		10x
DNA	10	—
10x Bsm	3	48
BSA	0.3	4.8
Bsm	0.2	3.2
H ₂ O	16.5	26.4
	30ul	20ul/tube

		4x
DNA	10	—
10x	3.0	1.2
BSA	0.3	1.2
Nco	0.2	0.8
Hind III	0.2	0.8
H ₂ O	16.3	6.52
	30.1	20ul/tube


3/19/96

Set up Digests of

1 HMSAF22 + G.P.PA2 (BglII/Asp)

2 HMSAF22 3' + G.P.PA2 (BglII/Asp)

	①	②
DNA	22	10
IOX	3	5
H ₂ O	221	23
BglII	1	1
Asp	1	1
	50ul	50ul

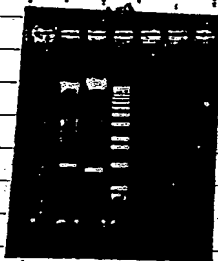Incubate 37°C
4 hrsRem. Load on Digest on gel with
11C ladderLooks like correct
digestion

Run on LMP Tomorrow

Make Probes: HTPAN08, HTEX07, HSH083
HFDJ02

	SAM	POS	CH	CPM	2SIG%	TIME
workup	HTPAN08	1	15	1	550513.31	0.49
	HFDJ02	2	16	1	945610.00	0.46
	HTEX07	3	17	1	761804.00	0.46
	HSH083	4	18	1	813240.00	0.50

3/20/96

screen HTO Lib

3 screens
all from CLF

Re Run Digests from 3/18 (pg 104)

3/20/96


1. HE2PM21 - Bam/Vho
2. HMSAF22 - Kpn/Xho
3. HSKBN01 - Kpn/Xho
4. HT45B02 - Bam/Vho
5. HTPAN08 - Nco/Htt

PC DNA has 2 Kpn
 Scaits to give 886 bp
 fragment

HT45B02 looks correct
 HSKBN01 looks correct

Need more Digests of HE2PM21 HTPAN08
 HMSAF22
 HT45B02

Run HMSAF22 on LMP Gel with
 1 Kb ladder.
 loaded samples after Carrie started Run.
 Cut out fragment
 Take picture

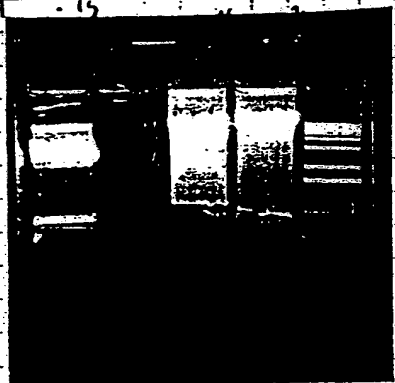

Gel slice ready
 for clean up
 or ligation


45H885

this
is CLF

24)

3/20/96

Run Digests of Mini preps.

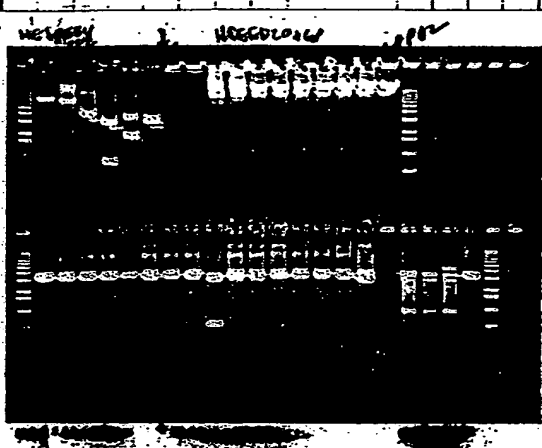

Bar Digests
 look like they
 are digested, but
 did not "pop out"
 fragment
 except HCEED20 #3.
 looks correct.

HCEED20 + PCEC
 None, look correct.
 No insert?

HCEED20.
 Should give ~400bp
 fragment that was
 digested w/ 600bp
 Ask J. J. J.

HCEED20 - Hard to tell
 Run on 0.8% or lower gel

Wash. HTO filters.
 0.2xSSC
 0.1% SDS. 3X. on film
 65°C, 1 hr.

- Lab meeting

01

5

60

60

65


10

1/2 hr


select
 in
 3/20

3/22/96

CORE NORTHERN BLOT SUBMISSION FORM

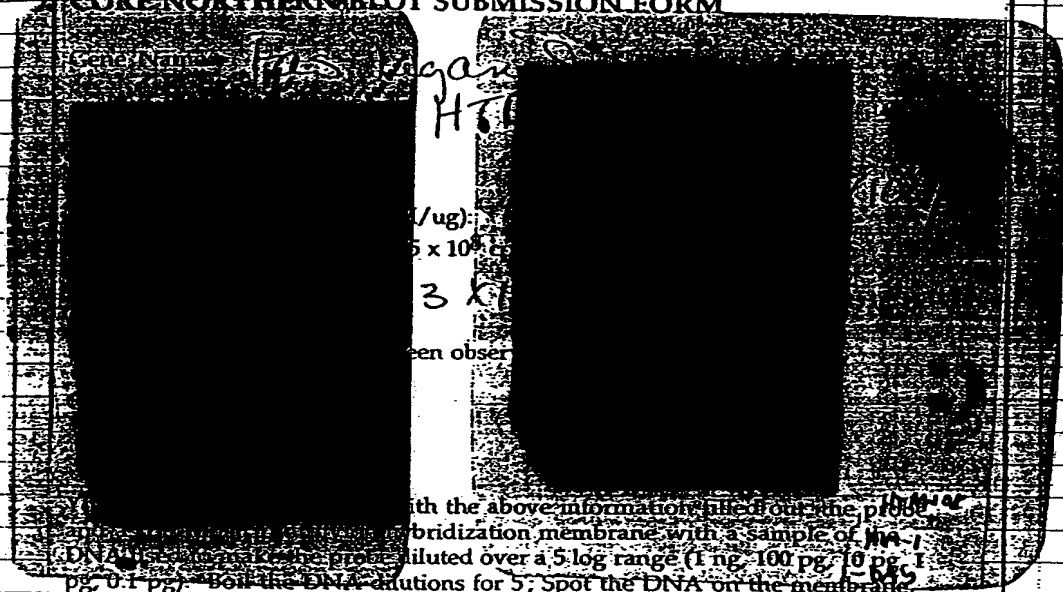

with the above information, filled out the probe hybridization membrane with a sample of DNA (1 ng, 100 pg, 10 pg, 1 pg, 0.1 pg). Boil the DNA dilutions for 5'. Spot the DNA on the membrane, using no more than 5 ul per spot. Alternately, denature the samples in Southern denaturation solution. In either case, crosslink the strip.

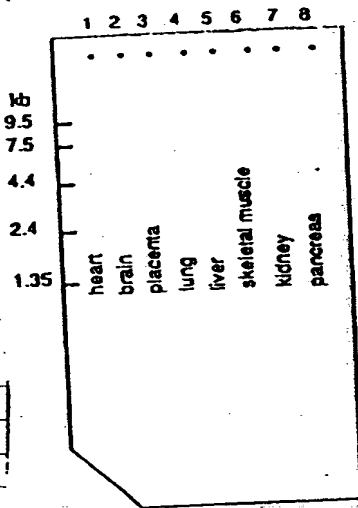

gram of your Northern band intensity of your band. We w

REA

5^{H5} , $I_n(5^{H5})$


ol

30.3×10^5

Ext. SHOWS NOTHING.

19-503, 15', 42°C



2x, 2x500, 1% 503, 30' EXCH, 65°C

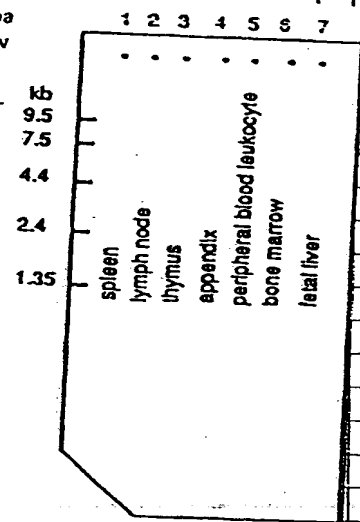

TNT Results

4/2/96

TNT RESULTS 4/2/96		
INVESTIGATOR	SAMPLE NAME	OBSERVED SIZE (KOs)
REINHARD EBER	HOECH07	51.8
YAJUN CHEN	HPMDH18 FRAME 1	48
YAJUN CHEN	HPMDH18 FRAME 2	48
YAJUN CHEN	HPMDH18 FRAME 3	48
ANN KIM	HTPAN08 + PCDNA	30
ANN KIM	HTPAN08 + 3HA	30, 35
ANN KIM	HT4SB02 + PCDNA	24
ANN KIM	HT4SB02 + 3HA	24
ANN KIM	HSBAW14 + GPPA2	24
CHARLES FLORENCE	HBMSE03	10
T7 POSITIVE CONTROL	DNASE 02-105	33
T7 NEGATIVE CONTROL	NO DNA	NONE
T3 POSITIVE CONTROL	HCAC183	33
T3 NEGATIVE CONTROL	NO DNA	NONE
REACTIONS PERFORMED BY: CARRIE FISCHER		

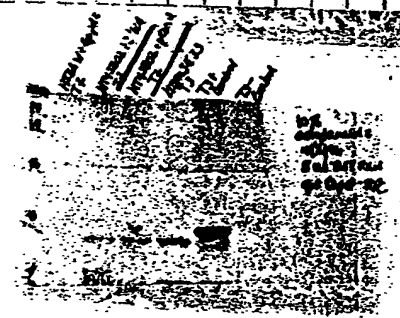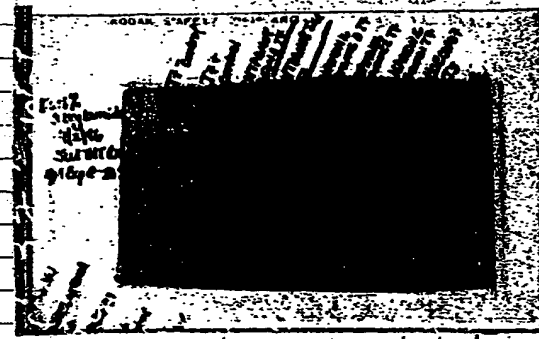

Made Probi for Northern MS.

HTPAN08
HTACB012
HTTB050
HTTB053

SAM	POS	CH	CPM	2SIG%	TIME
1	296	1	545253.31	0.49	0.30
2	297	1	930665.00	0.49	0.20
3	298	1	883680.00	0.48	0.20

CORE NORTHERN BLOT SUBMISSION FORM


149

4/14/96

Gene Name: Fas Ligand
 Clone Name: 1861-117RANO8
 cDNA size (in kb): [REDACTED]
 Specific Activity: [REDACTED]
 (MINIMUM ACTIVITY: [REDACTED] pg)
 Total counts: [REDACTED]
 Libraries in which: [REDACTED]
 Other: [REDACTED]

You must give us this form with the above information and a control strip of nylon hybridization membrane DNA used to make the probe diluted 1:10 (log 1 pg, 0.1 pg). Boil the DNA dilutions 5 min. Spot it using no more than 5 ul per spot. Alternatively, use Southern denaturation solution. In either case, c

You will receive the autoradiogram of your Northern blot 3 to 4 days, depending on intensity of your band. We will mark the film with the sizes of the marker ladder.

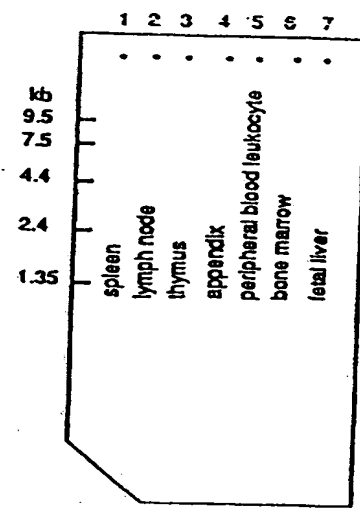

QC information

Name of hybridizer: Mark Porter
 Date hybridized: 4-3-96
 Blots hybridized: 1, #1
 Hybridization Solution: Hybrisol
 Counts/ml hyb. buffer added: 32.7 x 10⁵
 Exposure time: 1-DA
 Wash conditions: .2555c/.270 SDS

proax
4/5/96

1x, 15', 42°C
 2x, 20', 65°C
 PPT w/ EtOH
 pin lulanap
 Red OD 260/280

Sample ID	abs 260.0 nm	abs 280.0 nm	bkg abs 320.0 nm	abs 260.0 nm	abs 280.0 nm
1 H15AF22 168192	0.1731	0.1134	-0.0008	1.5122	
2 H15AF22 1680	0.1488	0.0968	-0.0026	1.5240	

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☒ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.