



Steven M. Ruben
Appl. No. 10/662,429

Department _____
Subject _____
Name Kim #7
Address 102894

 **National Brand** 43-648



Computation Notebook
Dennison Stationery Products Co., Framingham, MA 01701

 75 Sheets
11 1/4" x 9 1/4"
4x4 Quad.

0 73333 43648 8

BEST AVAILABLE COPY

Ruben EXHIBIT #90

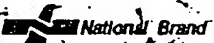Department _____
Subject _____
Name Kim #7
Address 102894
 43-648
Computation Notebook
Dennison Stationery Products Co., Framingham, MA 01701

75 Sheets
11" x 9"
4x4 Quad.
0 73333 43648 8

Ruben EXHIBIT 2090
Ruben v. Wiley et al.
Interference No. 105,077
RX 2090

TNT- HE9MF73/HTPAN/HTACI56/HPDDM9333

8PSM4

Do TNT on

(1) HE9MF73S05 0.25 μ g/ μ l
 (2) HE9MF73S07 0.25 μ g/ μ l
 (3) HTPAN08S04 0.25 μ g/ μ l
 (4) HPDDM93 0.25 μ g/ μ l
 (5) HTTCL56 0.1 μ g/ μ l

	(1)	(2)	(3)	(4)	(5)
Rabbit Lysozyme	25	25	25	25	25
Buffer	2	2	2	2	2
Pol T3	1	1	1	1	1
AA mix \ominus Met	1	1	1	1	1
³⁵ S meth.	1	1	1	1	1
RAasin	1	1	1	1	1
DNA	4	4	4	4	1
H ₂ O	15	15	15	15	10
	50	50	50	50	50

Incubate 30°C 1 1/2 hrs

5 μ l of TNT + 20 μ l 1X SDS Buffer
 Heat 95°C 2 min

Chill on ice

Quin Spin

Run 10 μ l on 12.5% SDS gel
 with Rainbow Marker. C+

Fix gel 20% Methanol 10% Acetic Acid
 30 min at 37°C

Amplify 30 min 37°C

Dry gel 70°C 2 1/2 hrs

1Bq/100

9/18
 03

TNT- HE9MF73/HTPAN/HTTACI56/HPODM9333

8/25/94

Do TNT on

EC
for and

(1) HE9MF73S05 0.25 µg/µl
 (2) HE9MF73S07 0.25 µg/µl
 (3) HTPAN08504 0.25 µg/µl
 (4) HPODM93 0.25 µg/µl
 (5) HTTACI56 0.1 µg/µl

Reagent	(1)	(2)	(3)	(4)	(5)
Lysate	25	25	25	25	25
Buffer	2	2	2	2	2
Pol. T3	1	1	1	1	1
RA mix + Met	1	1	1	1	1
35S meth.	1	1	1	1	1
RA prim	1	1	1	1	1
DNA	4	4	4	4	1
H ₂ O	15	15	15	15	10
	50	50	50	50	50

Incubate 30°C 1 1/2 hrs

5 µl of TNT + 20 µl 1X SDS Buffer
 Heat 95°C 2 min

Chill on ice

Quick spin

Run 10 µl on 12.5% SDS gel

with Rainbow Marker C+

Fix gel 20% Methanol 10% Acetic Acid
 30 min at 37°C

Amplify 30 min 37°C

Dry gel 70°C 2 1/2 hrs

ECO 1Bq/µl
976

see pg 03

TNT - HE9MF73/HTPAN/HTTACI56/HPDDM93

8/25/94

D. TNT on

TE.
1hr and

(1) HE9MF73S05 0.25 μ g/ μ l
 (2) HE9MF73S07 0.25 μ g/ μ l
 (3) HTPAN08S04 0.25 μ g/ μ l
 (4) HPDDM93 0.25 μ g/ μ l
 (5) HTTACI56 0.1 μ g/ μ l

p149
 book #108
 smk #6

Rabbit Lysozyme
 Buffer
 Pol T3
 AA mix @ Met
 35S meth.
 RNAse
 DNA
 H₂O

(1)	(2)	(3)	(4)	(5)
25	25	25	25	25
2	2	2	2	2
1	1	1	1	1
1	1	1	1	1
1	1	1	1	1
4	4	4	4	10
15	15	15	15	9
50	50	50	50	50

with

Incubate 30°C 1 1/2 hrs

Eco / Bgl II
976

5 μ l of TNT + 20 μ l 1X SDS Buffer
 Heat 95°C 2 min
 Chill on ice
 Quin Spin
 Run 10 μ l on 12.5% SDS gel
 with Rainbow Marker. C+

see pg
 03

Fix gel 20% Methanol 10% Acetic Acid
 30 min at 37°C
 Amplify 30 min 37°C
 Dry gel 70°C 2 1/2 hrs

34

TNT

8/25/94

Exposure:
-80°Cin Small
O/N

Cassette

8/26/94

Develop film

HE9MF73S03

HE9MF73S07

HTPAN08S04

HPDDM913

HTTC156


Expected Size

> ~ 83kd

~ 30.9 kd

~ 82.4 kd


~ 30.9 kd

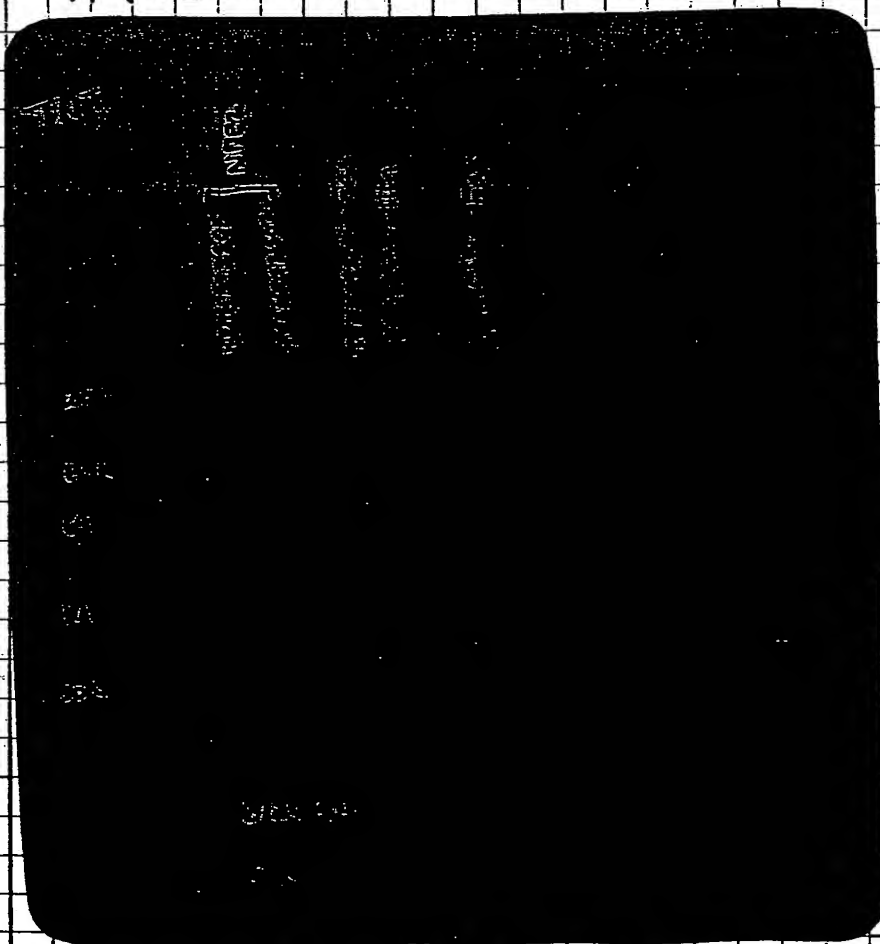HTPAN08
PO 140

Mini Pups / Maxi Pup

37

8/18/94

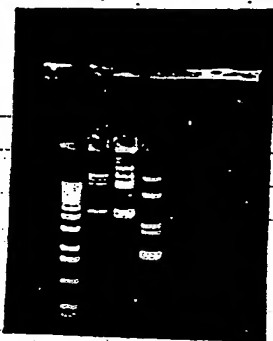

Sequence F & Rev.

Inoculate for maxis

HTTCIS6 B

HTPB411 C

HUNAD13 C

HCTBA03 A

HMSB243 A

HOSDK13 A

200 ml ~~ex~~
TBT Amp.

Incubate 37°C O/N

8/19/94

Obtain Maxi pup.

Spin Cells 4K 10min

Remove Supernatant

Resuspend pellet 10ml PI

Add 10ml P2 mix

Incubate RT for 10min 5min

Add 10 ml Chilled P3

Mix and incubate 20 min on ice

Spin 9K 30min

Equilibrate Tip 500 with 10ml QBT

Apply Supernatant through Kimwipe

Wash 2x 30ml QC

Elute 15 ml QF

Albert Rile
8/19/94

44

HE20142

(fragrant pup ^{HTPAN08504} HPRCC91)

8/31/94

Transfer Supernatant to fresh tube.

Run 1ul on gel with 1kb ladder.

HPRCC91 is wrong!
Try digesting with
Xho IHE20142 Xho I looks
good - ~600bp
use to probe filters

Digest	HPRCC91	Boiling minis - A+B
DNA	20ul	
10X #2	5	incubate 37°C
Eco RI	1	
Xho I	1	
H ₂ O	23	
	<hr/> 50ul	

5ul
pg 46
U

Digest	HTPAN08504	(Fas ligand full length)
DNA (0.6ug/ul)	15ul	
10X #2	10ul	
Eco RI	1ul	incubate 37°C
Xho I	1ul	
H ₂ O	73ul	
	<hr/> 100ul	

HE2 Continued
pg 47


46

Fragment prep HTPAN08S04 / HPRCC91

8/31/94

HTPAN08-P334

Run 2ul of digest on gel

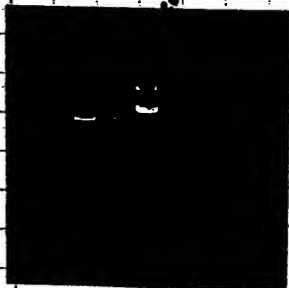

Does not look like
we have complete
Digestion.
Continue to digest
O/N at 37°C

9/1/94

Check Digests

Run 2ul on 1% gel with 1kb ladder.



HTPAN08S04 looks like it is
almost fully digested.

HPRCC91 does not look digested
at all.

Clean up DNA and try again.

2x phenol extract w/ equal volume
2x chloroform extract w/ equal volume
Ethanol / 3M NaAcetate precipitate DNA
1x 70% Ethanol wash
Dry pellet
Resuspend in 50ul TE

Set up

Digest

DNA

20ul

10x#2

5ul

H₂O

23ul

Xho

1ul

EcoRI

1ul

50ul

Incubate at 37°C.

P334

HTPAN08504 - Protein - fragment Prep. 157

9/2/94

PCR -

HTPAN08504 1 μ g μ l.

7690 3' Asp. 718	0.3	12x
7689 5' BamHI	0.3	3.6
10x PCR	3.5	3.6
10x dNTP	3.5	42
Taq	0.2	42
H ₂ O	26.2	2.4
DNA	1	314
	35 μ l.	12
		420

Should
get a
fragment
~ 850 bp.

PCR Program # 69

95°C 5min
98°C 20 sec
55°C 20 sec
72°C 1min
72°C 7 1/2 min
4°C Hold

30x

Run 5 μ l on gel

looks good -

Add equal Volume of
13% PEG / 1.6M NaCl
let sit on ice 10 min

Spin 10 min
Remove Supernatant
Add 1ml 70% Ethanol
Spin 5 min
Remove Supernatant
allow pellet to dry
Resuspend in 100 μ l TE
Store -20°C

HTPANO8504 - Fragment Prep

9/6/94

Digest fragment - $\text{Hsp 718} + \text{Bam HI}$
 $\text{Hsp Buffer B}, \text{Bam} - \text{Bam Buffer}$

DNA (PCR product)	20 μl
10X Buffer	5 μl
H ₂ O	24 μl
Enzyme	1 μl
	<hr/> 50 μl

Incubate 37°C 4 hrs.

Add Complement Enzyme.
 Incubate 37°C 4 hrs.

Run on 0.8% LMP Agarose gel
 with 1 kb ladder.



Cut out fragments into
 1.5 ml Tube

Add 400 μl NaI
 Heat 55°C 5 min

Mix

Add 5 μl Glass milk

Vortex

Let sit at RT 5 min

Spin 5 sec

Remove supernatant

Resuspend pellet in 200 μl
 & Wash Solution

3x Spin 5 sec

Remove Supernatant

Spin 5 sec

Remove Supernatant

Add 20 μl TE

Heat 55°C

Spin 5 sec

1 min


(pg 69)

68

FRAGMENT Prep HPRCC91 / HIPANOS8504 Asp / Bam.

(pg66) (pg 69)

9/7/04

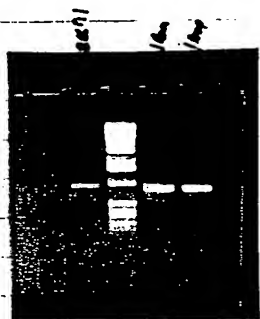Run 1 μ l on gel with 1 kb ladderHPRCC91 Kpn / Bam HI
 ≈ 100 ng/ μ lHIPANOS8504 3' Asp / 5' Bam
Asp / Bam Digest
Gene Cleaned 9/6
 ≈ 300 ng/ μ lHIPANOS8504 3' Asp / 5' Bam
Bam / Asp Digest
Gene Cleaned 9/6
 ≈ 300 ng/ μ lStore in Fragment Box -20°C .

HTPANC8504 - PAZ Protein Prep

69

(pg 58)

9/6/94

Transfer Supernatant to fresh Tube.
 Resuspend pellet in 10ul TE
 Heat 55°C 1 min
 Spin 5.8cc

Transfer Supernatant to fresh tube.

Set up ligations w/ PAZ Asp/Asn digest - (See pg 67)

①	
HTPANC8504 Asp/Asn	3
PAZ Asp/Asn	1
10x Buffer	1
H ₂ O	14
T4 ligase	1
	20ul

②	
HTPANC8504 A/B	3
PAZ B/A	1
10x Buffer	1
H ₂ O	14
T4 ligase	1
	20

③	
HTPANC8504 B/A	3
PAZ A/B	1
10x Buffer	1
H ₂ O	14
T4 ligase	1
	20

④	
HTPANC8504 B/A	3
PAZ A/A	1
10x Buffer	1
H ₂ O	14
T4 ligase	1
	20

⑤	
HTPANC8504 A/B	3
PAZ B/B	1
10x Buffer	1
H ₂ O	15
T4 ligase	1
	20

⑥	
HTPANC8504 B/A	3
PAZ B/B	1
10x Buffer	1
H ₂ O	15
T4 ligase	1
	20

⑦	
Fragment PAZ A/B	1
10x Buffer	1
H ₂ O	17
T4 ligase	1
	20

⑧	
PAZ B/A	1
10x Buffer	1
H ₂ O	17
T4 ligase	1
	20

⑨	
Fragment PAZ B/B	1
10x Buffer	1
H ₂ O	18
T4 ligase	1
	20ul

Incubate 16°C O/N

(See pg 68 for fragment conc)

70

HTPANO8504 Problem. pA2

9/7/94

Thaw DH5 α Chemically Competent Cells
on ice.

Aliquot 100 μ l of Cells into fresh tubes
Add 10 μ l of Restriction Rxd.

Mix gently
Incubate on ice 1 hr.

Heat 42 $^{\circ}$ C 45 sec

place on ice

Add 400 μ l LB

incubate 37 $^{\circ}$ C 1 hr.

Plate onto LB + Amp Plates.

- ① HTPANO8504 Asp/Bam + pA2 Asp/Bam - 50 + 200 μ l
- ② HTPANO8504 Asp/Bam + pA2 Bam/Asp - 50 + 200 μ l
- ③ HTPANO8504 Bam/Asp + pA2 Asp/Bam - 100 μ l
- ④ HTPANO8504 Bam/Asp + pA2 Bam/Asp - 100 μ l
- ⑤ HTPANO8504 Asp/Bam - 200 μ l
- ⑥ HTPANO8504 Bam/Asp - 200 μ l
- ⑦ pA2 Asp/Bam - 200 μ l
- ⑧ pA2 Bam/Asp - 200 μ l
- ⑨ Lig Run only - 200 μ l
- ⑩ 5ng pA2 plasmid DNA - 200 μ l
- ⑪ DH5 α Cells only - 200 μ l

incubate 37 $^{\circ}$ C O/N.

9/8/94

Inoculate Colonies into LB + Amp in

96 Well dish
Incubate at 37 $^{\circ}$ C w/ aeration 4 hrs

PCR to check inserts

HTPAN08 SOY

Protein

PAZ

71

PCR



① #7689	1	(22x)	22
#7690	1		22
10x dNTP	3.2		70.4
10x PCR	3.2		70.4
H ₂ O	21.4		470.8
Taq	0.2		4.4
Culture	2		—

Total 32 μ l

PCR Program #69.

95°C 5min
 95°C 20sec
 55°C 20sec
 72°C 1min
 72°C 7 1/2 min
 4°C Hold.

30x

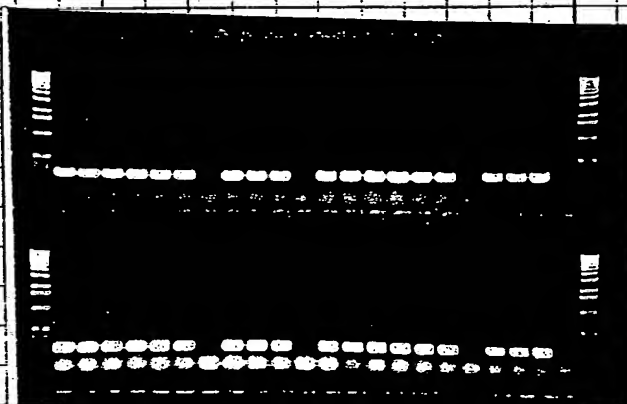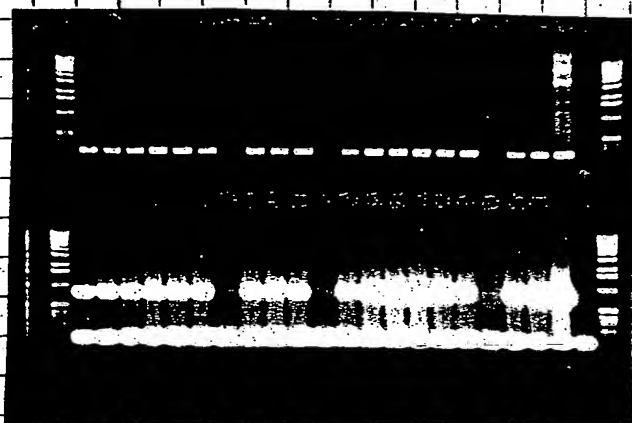Run 10 μ l on gel with 1 Kb ladder.

There are clearly positives.
 Do boiling mini + Sequence with internal primers
 out through ligation junction

72 HTPAN08504 into pA2.

9/8/94

Inoculate 5ml TB + Amp with all the
+ clones -

Incubate 37°C O/N w/aeration

9/9/94

Boiling mini preps

Spin 2ml of culture in 2ml tubes
6 min

Remove supernatant

Resuspend pellet in 450 μ l of
STE1 + RNase + Lysozyme

Boil 1 min

Spin 10 min

Remove pellet w/ toothpick

Add 75 μ l of 13% PEG-8000 / 1.6M NaCl

Vortex well

Spin 10 min

Remove supernatant

Add 1000 μ l 70% Ethanol

Mix

Spin 5 min

Remove supernatant

Allow pellets to dry

Resuspend in 150 μ l of TE

Run 2 μ l on gel with 1 kb ladder
and pA2 (control)

BD1 \rightarrow BD 17

Add 150 μ l more TE

Digest to make sure
insert "pA2" out



pg 81

HTPAN08S04 - pA2

Boce. Virus.

81

(p72)

9/17/99

HTPAN08B01 - 17


Digest with Pstm/Asp.

DNA	5 μ l	(17x)
10x B	3	51
H ₂ O	21.6	367.2
Bam/Asp	0.2 / 0.2	3.4 / 3.4
	30 μ l.	

Incubate 37°C 2 hrs

Run 10 μ l on gel with 1kb ladder

They all look correct.

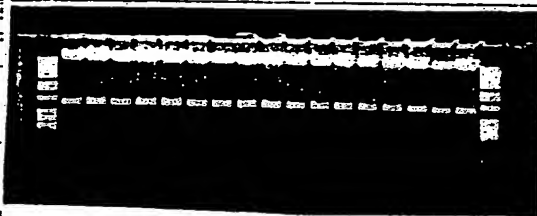

Submit for
Sequencing
to confirm
Sequence and
to confirm Cloning
sites are OK

Submit in 96 Well dish

HTPAN08S04

HTPAN08B01

B02

B03

B04

B05

Use Primers:

RP01

EP03

RP05

RP06

EP07

RP12

EP14

EP16

FA7

FA8

freeze Glycerols of B01-B17 -

9/12/94

Incubate 150 ml TBT Amp
w/ glycerol stock of
HTPAN08801
HTPAN08802

Incubate 37°C O/N w/aeration

9/13/94

Quagen Maxi Prep

Spin Cells 4500 rpm 15 min
pour off supernatant
Resuspend pellet in 10 ml P1 + RNase
Add 10 ml P2 mix well
let sit RT 5 min
Add 10 ml cold P3 mix
well

let sit on ice 20 min

Spin 4500 rpm 15 min

Equilibrate Tip 500 with 10 ml DBT
Apply Supernatant through Kimwipe
Wash Column 2x
with 30 ml QC Buffer

Elute DMT 15 ml QF

Add isopropanol 0.7 volumes (10.5 ml)

Mix Well

Spin 9K GSA Rotor 30 min

Pour off Supernatant

Add 15 ml 70% Ethanol (ice
cold)

Spin 15 min 9K GSA Rotor

i Prep

HTPAN08504-PAZ Maxi Prep.

88

9/13/94

Carefully Pour off Supernatant
 Allow pellet to air dry at
 Room Temp ~ 30 min

Resuspend pellet in 400 μ l of TE
 Transfer to 2 ml Tube
 Read OD_{260/280}: 2.00
 Dilute
 DNA 5 μ l \rightarrow 99.5 μ l H₂O

Sample	abs		abs	avg abs	260.0 nm		280.0 nm
	260.0 nm	280.0 nm			260.0 nm	280.0 nm	
1201	0.1179	0.0760	0.0188	1.7324	0.5772	1.16 μ g/ μ l	472 μ g
2202	0.1309	0.0847	0.0170	1.6834	0.5940	1.31 μ g/ μ l	524 μ g

Store 4°C in Plasmid Box.
 Will need to do IP - Transfection

Run 0.5 μ l of plasmid + PAZ and
 Lambda Marker



too much λ marker, but
 clones look good

9/22/94

Give Jim 2
 20 μ l each plasmid for him
 to transfect Cos Cells.

as per Steve Ruben's Request (9/22/94)

(PA III)

HE20142

97

Plate ~~5~~ ⁵ on LB + Amp + IPTG / Xgal 9/27/94
 plate

9/27/94

Plates did not grow
 Try PCRing for issues.

① M13R + FP10 ~500 bp
 ② RPO1 + FPO9 ~430 bp
 ③ M13R + FPO8 ~1040 bp

		27x
M13R	0.1 μ l	3.7
FP10	1.3	35.1
10x dNTP	3.5	94.5
10x PCR	3.5	94.5
Taq	0.2	5.4
H ₂ O	21.4	577.6
Rescue	5	

		27x
M13R	0.1	2.7
FPO8	1.1	29.7
	3.5	94.5
	3.5	94.5
	0.2	5.4
	21.6	583.2
	5	

		27x
RPO1	0.2	5.4
FPO9	1.3	35.1
	3.5	94.5
	3.5	94.5
	0.2	5.4
	21.3	575.1
	5	

⊕ Control HE20142 plasmid
 ⊖ Control P20

Run PCR prog #69

90°C 5min
 95°C 20sec
 55°C 20sec } 32x
 72°C 1min
 72°C 7 1/2 min
 4°C Hold

Run Pcr on gel as - forgot to add
 run 1 kb ladder

Abner M. R. L.
 9/27/94

HTPB411S15 & HTPAN08504 - into PD10

111

pg 83

9/28/94

Fragment prep for HTPB411S15 3'x12" / 5'8"m
HTPAN08504 3'x12" / 5'8"m

DNA	HTPB411S15
10x dNTP	10
10x PCR	10
2500 2501	0.3
2502	0.3
Taq	0.3
H ₂ O.	78.1
	100ul.

DNA	HTPAN08504
10x dNTP	10
10x PCR	10
2499	0.3
2500	0.3
Taq	0.3
H ₂ O	78.1
	100ul.

Do 10 Runs of each and 1 negative control.

PCR Prog # 69
95°C 5 min
95°C 20 sec
55°C 20 sec } 30x
72°C 1 min
72°C 7 1/2 min.
4°C hold.

Run 5ul on gel with 1 kb ladder



- HTPAN08504 ~ 800 bp


- HTPB411S15 ~ 2400 bp

- PEG precipitate - equal Volume
Spin 10 min
Remove supernatant
Add 1 ml 70% ethanol
Spin 5 min
Remove Supernatant

Allow pellet to dry
Resuspend in 150ul of TE

9/28/94

Run on 0.8% LMP gel

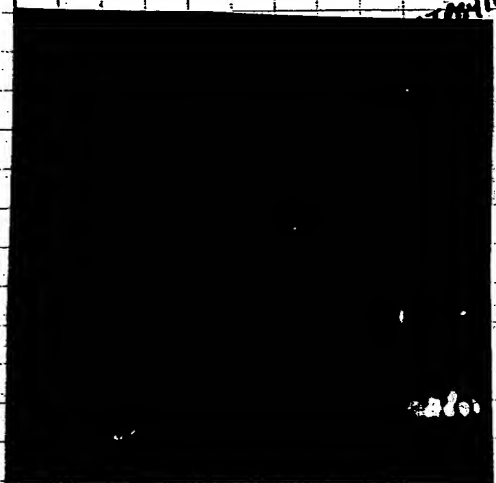Cut out of gel
Gene clean

Add 1ml NaI

Heat 55°C 5min

Mix

Add 7ml Glass Milk

let sit at RT 5min

with occasional
mixing

Spin 5 Sec

Remove supernatant

Resuspend pellet in
400ul Wash Buffer

3X

Spin 5 Sec

Remove supernatant

Spin 5 Sec

Remove supernatant

Add 30ul TE

Heat 55°C 2min

Spin 5 Sec

Transfer to fresh tubes

Add 20ul TE

Heat 55°C 2min

Spin 5 Sec

Transfer to tube

- Run 2ul on gel with 1 Kb ladder.

9/28/94

looks good
Digest with

Bam

DNA 30ul

10x#2 5

H₂O 13

Enzyme 1

5ul

Digest with Bam
for 4 hrs at 37°C
then add 1ul
KpnI & digest at
37°C for 4 hrs

store 4°C

HTPB411515 + HTPAN08504 in PD10

113

9/22/94

Run samples on 0.8% LMP gel



HTPAN08504 Barn/1X6
Digest
Cut out band 0.9 kb
And 1.6 kb fragments
HTPB411515 Barn/1X6
Digest
Cut out 2.4 kb frag

Gene Clean
Add 1 ml NaI
Heat 55°C 5 min
Add 7 ul Glass milk
mix & let sit at RT 2 min

Spin 5 sec
Remove supernatant
Resuspend pellet 400 ul Wash Buffer
Spin 5 sec
Remove supernatant
Spin 5 sec
Remove supernatant
Resuspend pellet 20 ul TE
Heat 55°C 2 min
Spin 5 sec
Transfer to fresh tube
Add Resuspend pellet in 20 ul TE
Heat 55°C 2 min
Spin 5 sec
Transfer to tube
Run 2 ul on gel with 1 kb ladder

1 - 0.9 kb HTPAN08504
2 - 1.6 kb HTPAN08504
3 - 2.4 kb HTPB411515



114

HTPAN08504 & HTPB411S15 in PD10

9/30/94

Set up ligations

- ① HTPAN08504 0.9 Kb Bam/Xba
 ② HTPAN08504 1.6 Kb Bam/Xba
 ③ HTPB411S15 2.4 Kb Bam/Xba

DNA Fragment + Vector

#1	①	4ul +	pD10 Bam/Xba	2ul
#2	①	4ul +	pD10 Xba/Bam	2ul
#3	②	6ul +	pD10 Bam/Xba	2ul
#4	②	6ul +	pD10 Xba/Bam	2ul
#5	③	4ul +	pD10 Bam/Xba	2ul
#6	③	4ul +	pD10 Xba/Bam	2ul
#7	①	4ul	+	—
#8	②	6ul	+	—
#9	③	4ul	+	—
#10	—	—	+	pD10 Bam/Xba
#11	—	—	+	pD10 Xba/Bam
#12	—	—	+	pBSK rest digest
#13	—	—	+	—

2ul Total Reaction Volume

① Control m5 cells only ② control 10ng pD10

Incubate RT 1 hr

Add 10ul to 100ul M15 cells

Let sit on ice 1 hr

Heat 42°C 45 sec

Sit on ice

Add 400ul LB

Let incubate 37°C 1 hr

plate 200ul on LB + Amp + Kan plates

for pD10 vectors

plate 100ul on LB + Amp + Kan plates

for pBSK

Let sit RT over the weekend.

HTPAW08504 + HTPB411515 in PD10

115

10/3/94

Picked ⑥

Clones of HTPA W08504 +
HTPB411515

into 200 μ l LB + Amp + Kan.
Do PCR. Incubate 37°C, 4 hrs.

HTPA W08504

HTPB411515

10x PCR	3.2
10x dNTP	3.2
2499	0.3
2500	0.3
100	0.2
H ₂ O	22.8
Culture	2 μ l
	32 μ l



10x PCR	3.2
10x dNTP	3.2
2501	0.3
2502	0.3
100	0.2
H ₂ O	22.8
Cult	2
	32 μ l

PCR Program # 69

95°C	5 min
95°C	20 sec
55°C	20 sec
72°C	1 min
72°C	7 1/2 min
4°C	Hold

④ Controls -
use 10 μ g Plasmid
DNA

Run 10 μ l on gel with 1 kb ladder

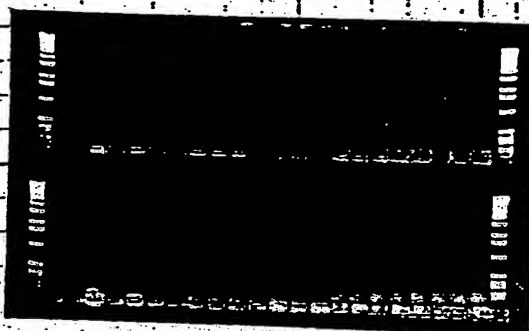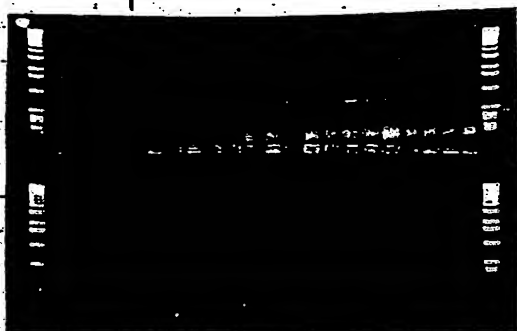

HTPA W08

10/3/94



HTPANOS



HTPB411

Transform Remaining Amount of Ligations
into M15 cells

Add DNA to 100 μ l M15 cells.

Incubate 1 hr.

Heat 42°C 45 sec

Place on ice

Add 400 μ l LB

Incubate 1 hr.

Plate onto LB + Amp + Kan Plates

Incubate O/N at 37°C

10/4/94

Pick Colonies ~~on~~ into LB + Amp + Kan

media in 96 well dish.

Incubate 37°C w/ aeration

Do PCR

See Run map from 115

See PG 123

HTRAXS - HTPBU

in PD10



123

P116


1/4/94

PCR Prog H69

Run 100 on Gel w/ 1 kb ladder

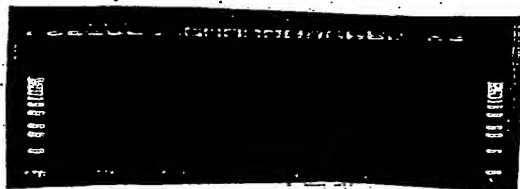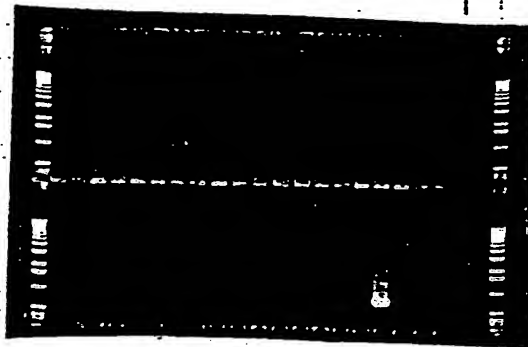



HTRAXS

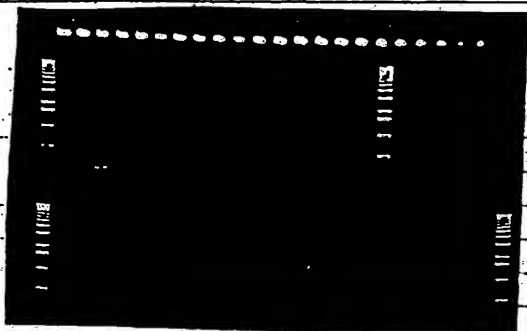

124

HTPAN08 or HTPB411 in pD10

10/4/94

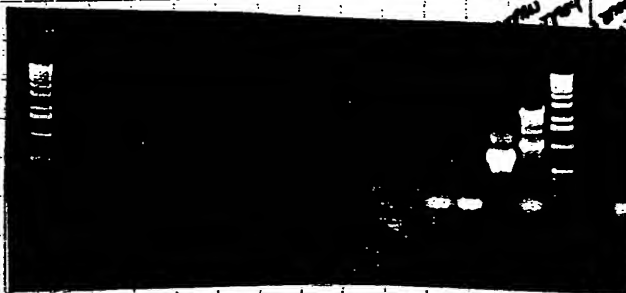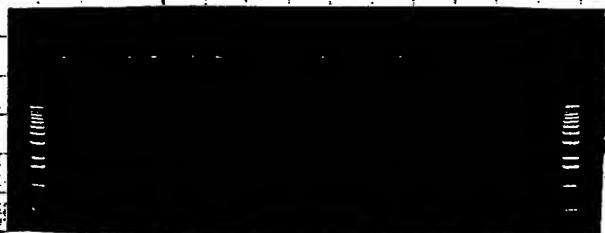

looks like 1 ϕ for HTPAN08
 Remata 1000 P10 Bam/Xba from Vector.
 Grow up Maxi Prep.

0/5/94

Re set up ligations

	1	2	3	4	5	
HTPAN08 0.9	6	—	—	—	—	✓
HTPAN08 1.6	—	8	—	—	—	✓
HTPB411 S15	—	—	8	—	—	✓
pD10 Bam/Xba	2	2	2	2	—	✓
10X	2	2	2	2	2	✓
14 ligase	1	1	1	1	1	✓
H ₂ O	9	7	7	15	17	

Incubate 16°C O/N.

0/6/94

Transform MIS cells

H1111111111 + P1111111111 in PD10

125

10/6/94

to 100ul M15 Chemically Competent
Cells add 10 ul of ligation
Incubate all on ice 1 hr
Heat 42°C 45 sec
Place on ice
Add 400ul LB
Incubate 37°C 1 hr
Plate 200ul into LB + Amp/Kan 150mm
plates
Incubate 37°C O/N

10/7/94

~~Observation~~
Colonies did not grow well
only about 10 colonies per plate.
Remake PD10 Vector.
Remake insert fragments

10/10/94

Digest PD10 1.57 ug/ul

DNA	3.2
10X	10
H ₂ O	84.8
Bam/Xba	1/1
	100ul

Incubate 37°C O/N

10/11/94

Run on 0.8% LMP Gel with 1Kb ladder
80V 1 1/2 hrs.



Isolate gel fragment into
Tube - 2

Add 1000 μ l NaI.
Heat 85°C 5min until
melted

Mix well
Add 7 μ l Glass milk
Mix well

Incubate at RT 5min
with occasional mixing

Spin 7sec

Remove supernatant

Resuspend pellet in 400 μ l
Wash Buffer.

3X

Spin 7sec

Remove Supernatant

Spin 7sec

Remove Supernatant

Resuspend pellet 20 μ l TE

Heat 85°C 1min

2X

Spin 10sec

Transfer & Combine into 1 fresh tube

Run 2 μ l on gel with 1Kb ladder.



Make more fragment
See pg 146

Hr2CC91 HPRDM93

137

10/13/94

Develop film

Pick ⊕ clones into 500ul Sp Buffer.
only HPRCC91 seemed to work -
115 ⊕ clones
HPRDM93 - Nothing lit up
Try HPRDM Again (See pg 147)

10/18/94

Do PCR on ⊕ clones from 1^o phage.

10x PCR	3.5	120x
10x dNTP	3.5	420
M13R	0.1	420
(8235) FRO3	2	12
1ug	0.2	240
H ₂ O	20.7	24
phage	5	2484
	35ul	30ul/tube

Run PCR Program # 69.

95°C 5min
95°C 20 sec
55°C 20 sec
72°C 1min
72°C 7 1/2 min
4°C Hold

30x

Re Use HPRCC93
plasmid DNA
as ⊕ control
Use H₂O as
⊖ control

Run 10ul on 1% Agarose gel with
1kb ladder

Neer M Ph.D.
10/15/94

146

HTPB411 + HTDAS05 into PD10

pg 126

10/11/94

Take PCR'd product from.
(9/23/94 pg 112)

Digest

DNA	3
10x Buffer	5
H ₂ O	49.0
XbaI/PstI	0.5/0.5
	<hr/> 50ul

Incubate 37°C 4 hrs

Run on 0.8% TAP gel

Excise Band

Gene Clean (See pg 126 for protocol)

Run 2x on gel



10/12/94

Set up digestion again.

With PD10 from 10/11 &
 HTDAS08 band from 10/11
 HTPB411 band from 10/11

HIPB411 & H7 PANO 8 with PD10 147

10/12/94

~~Exceeded 100 cells (H7)~~

	1	2	3	4	5	6	7	8
HIPANO 80.9	4			4				
HIPANO 81.6		4			4			
HIPB411			4			4		
PD10 B/xba	1.5	1.5	1.5			1.5	1.5	
10x T4 Buffer	2	2	2	2	2	2	2	2
H ₂ O	11.5	11.5	11.5	13	13	13	15.5	17
T4 ligase	1	1	1	1	1	1	1	1
	20	20	20	20	20	20	20	20

Incubate 37°C O/N.

10/13/94

Transform M15 cells.

Thaw M15 Chem. Competent cells on ice

To sterile tubes combine

100 µl M15 cells

10 µl of ligation

let sit at RT 1 hr

Heat 42°C 45 sec

Quick Chill

Add 400 µl LB

Incubate 37°C 1 hr

plate 200 µl into LB + Amp +

Kan plates (150 mm)

Incubate 37°C D/N.

use 10 µg PD10
as @ Control.

HTPB411 & HTPAN08 into PD10

10/14/94

Plates look OK -

Plates with PD10 Vector alone
has many colonies

Pick 100 of each clone into

200 μ l of LB + Amp^r Kan

Incubate at 37°C w/ aeration

PCR.

HTPAN08	
2499	0.2
2500	0.2
10x dNTP	3.2
10x PCR	3.2
Taq	0.2
H ₂ O	2.3
Culture	2
	<hr/> 32 μ l

HTPB411	
2501	0.2
2502	0.2
10x PCR	3.2
10x dNTP	3.2
Taq	0.2
H ₂ O	2.3
Culture	2
	<hr/> 32 μ l

Run Program #69.

95°C	5 min	} 30x
95°C	20 sec	
55°C	20 sec	
72°C	1 min	
72°C	7 1/2 min	

Run 10 μ l on gel with 1 kb ladder -Nothing Showed up!

Try Again

H1PB411 & HTPANOS into PD/O

249

10/19/94

PCR Amplify fragments Again

① HTPANOS

2499	0.2
2500	0.2
10x PCR	10.0
10x dNTP	10.0
Taq	0.3
ONAP buffer	1 (100µl)
H ₂ O	78.3
	<hr/> 100µl

② H1PB411

2501	0.2
2502	0.2
10x PCR	10
10x dNTP	10
H ₂ O	78.3
Taq	0.2
ONAP	1
	<hr/> 100µl

Set up 5 tubes of each
Run PCR

95°C	5min	30x
95°C	20sec	
55°C	20sec	
72°C	1min	
72°C	7 1/2 min	

Run 10 µl on gel with 1 kb ladder



Precipitate PCR Prep
Add equal volume
13% PEG / 1.6M NaCl
Let sit on ice 15 min
Spin 15 min

1x 70% ethanol Wash
Dry pellet slightly
Resuspend in 50 µl TE

150

H18B11 & H18A08 into pD10

10/19/94

Run gel on 1% TAE gel with 1 Kb ladder


looks good -
Set up Digestions

DNA	40ul
H ₂ O	48ul
10x #2	10
Boa/Boa	1/1
	<hr/> 100ul

Incubate 37°C O/N

10/20/94

Run on 0.8% LMP Gel with 1 Kb ladder

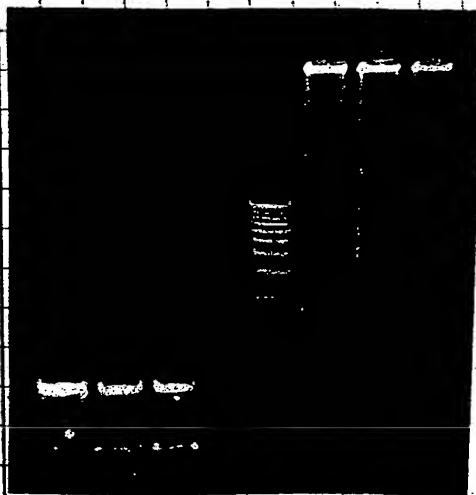Gene Clean
Seq Page 120

Boa

Set up sequencing

with pD10 from
10/19 - Seq 151

Seq 3 Boe 3

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☐ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☒ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.